


# Immunology Research Review

An Overview of Recent Immunology Research Publications Featuring Illumina® Technology

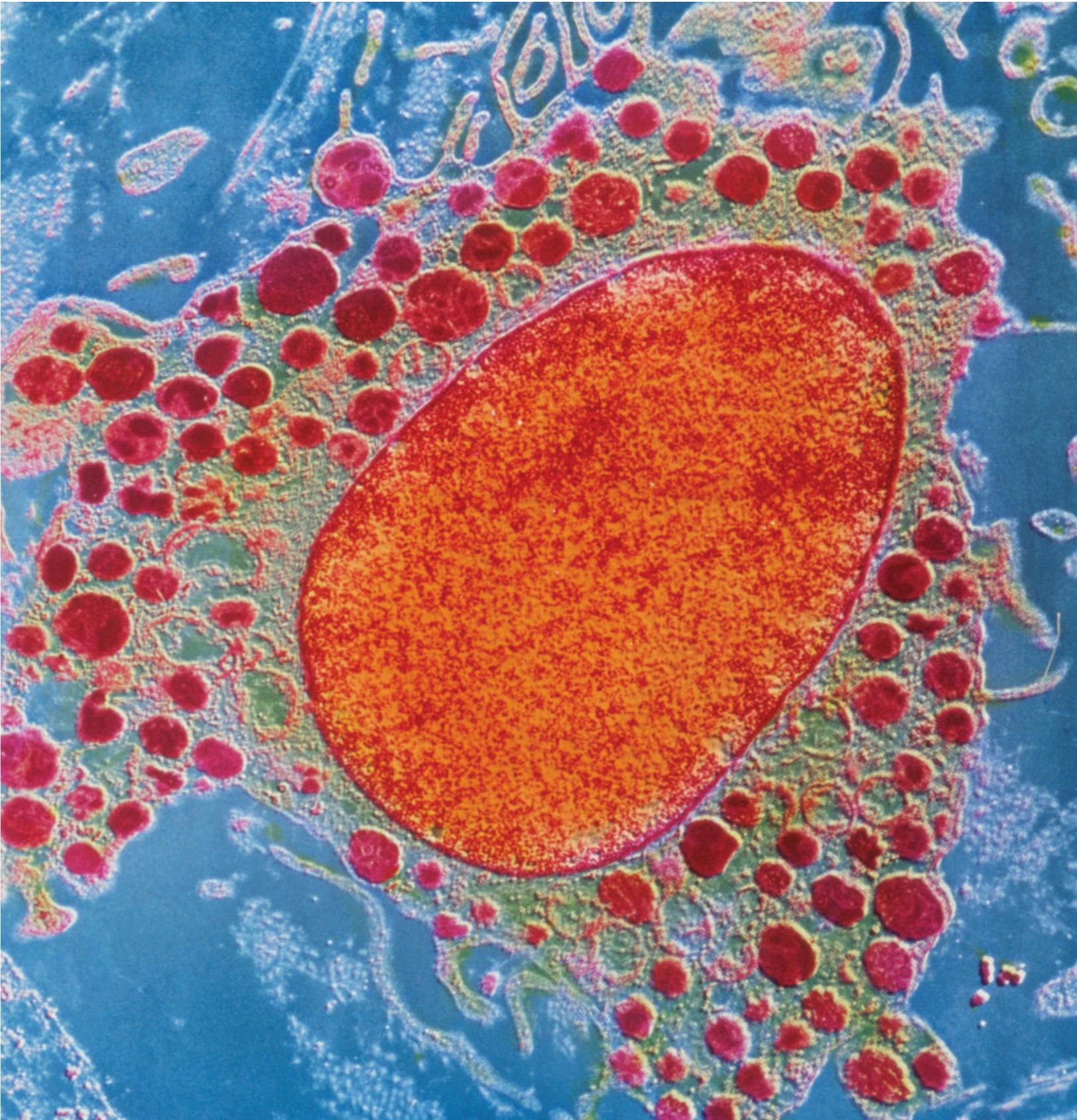

# CONTENTS

- 4** Introduction
- 6** Adaptive immunity
  - Repertoire Sequencing of Lymphocyte Receptors
  - Single-Cell Repertoire Sequencing
- 13** Lymphocyte Development
  - T Cell Development
  - B Cell Development
- 18** Innate Immunity
- 20** Cancer and the Immune Response
- 21** Microbiota and the Immune System
- 24** Major Histocompatibility Complex
  - Phase-Defined HLA Sequencing
- 29** Self vs Non-Self Antigen Discrimination
  - Tolerance
  - Autoimmunity
  - Solid Organ Transplantation
- 36** Infectious Diseases and Vaccines
  - Viral Infections
  - Vaccine Development
- 41** Techniques
  - miRNA and noncoding RNAs
  - ChIP-Seq
- 46** Bibliography

This document highlights recent publications that demonstrate the use of Illumina technologies in immunology research. To learn more about the platforms and assays cited, visit [www.illumina.com](http://www.illumina.com).

# INTRODUCTION


Immunology is the field of study concerned with the recognition and disposal of foreign or “non-self” material that enters the body. This material is usually in the form of life-threatening infectious microorganisms<sup>1</sup> or cancer<sup>2</sup> but sometimes, unfortunately, in the shape of life-saving graft transplantation.<sup>3</sup> The body can also be tricked into mobilizing the immune response against itself, to create autoimmune diseases. The NIH estimates that approximately 23.5 million Americans suffer from autoimmune disease and that the prevalence is rising.<sup>4</sup> Recent progress in the treatment of autoimmune diseases amply illustrates the impact that advancements in immunology are having on human health and disease.

Next-generation sequencing technology is proving to be a powerful tool to map the vast repertoire of immune cells that are capable of recognizing the seemingly boundless array of targets.<sup>5</sup> Repertoire sequencing has enabled researchers to identify unique receptor variants found in individuals with susceptibility to hematological malignancies, autoimmune diseases, and allergen response.<sup>6</sup> This approach is rapidly gaining the attention of translational scientists who seek to improve patient care. Hematologists have led the repertoire sequencing effort and have demonstrated the reliability, cost-effectiveness, and medical value of repertoire sequencing in hematopoietic stem cell transplantation.<sup>7</sup>

The major histocompatibility complex (MHC) is a locus that encodes a highly variable repertoire of cell surface proteins that present foreign antigens to T-cells. The encoded repertory of cell-surface molecules enables immune recognition and clearance of foreign agents. Genes within this locus are routinely assessed in matching patients and donors for solid organ transplantation<sup>8</sup> and hematopoietic stem cell transplantation.<sup>9</sup> By comparing variants of these genes between healthy and affected individuals, researchers are now able to elucidate the root causes of disease susceptibility (i.e. hematological, autoimmune, allergies, hypersensitivities, chronic inflammation, infectious diseases).<sup>10</sup>

The development of Illumina’s next generation sequencing provides the quality, throughput and read lengths required by the research community to map the human immune response at high resolution. The emergence new approaches such as phase-defined sequencing and single-cell sequencing can be expected to accelerate this knowledge base.

- 
1. Neller M. A., Burrows J. M., Rist M. J., Miles J. J. and Burrows S. R. (2013) High frequency of herpesvirus-specific clonotypes in the human T cell repertoire can remain stable over decades with minimal turnover. *J Virol* 87: 697-700
  2. Haen S. P. and Rammensee H. G. (2013) The repertoire of human tumor-associated epitopes--identification and selection of antigens and their application in clinical trials. *Curr Opin Immunol* 25: 277-283
  3. Meyer E. H., Hsu A. R., Liliental J., Lohr A., Florek M., et al. (2013) A distinct evolution of the T-cell repertoire categorizes treatment refractory gastrointestinal acute graft-versus-host disease. *Blood* 121: 4955-4962
  4. The Cost Burden of Autoimmune Disease: The Latest Front in the War on Healthcare Spending. 1–14 ([www.diabetesed.net/page/\\_files/autoimmune-diseases.pdf](http://www.diabetesed.net/page/_files/autoimmune-diseases.pdf))
  5. Woodsworth D. J., Castellarin M. and Holt R. A. (2013) Sequence analysis of T-cell repertoires in health and disease. *Genome Med* 5: 98
  6. Robins H. (2013) Immunosequencing: applications of immune repertoire deep sequencing. *Curr Opin Immunol* 25: 646-652
  7. Warren E. H., Matsen F. A. t. and Chou J. (2013) High-throughput sequencing of B- and T-lymphocyte antigen receptors in hematology. *Blood* 122: 19-22
  8. Eapen M., Rubinstein P., Zhang M. J., Stevens C., Kurtzberg J., et al. (2007) Outcomes of transplantation of unrelated donor umbilical cord blood and bone marrow in children with acute leukaemia: a comparison study. *Lancet* 369: 1947-1954
  9. Marks C. (1983) Immunobiological determinants in organ transplantation. *Ann R Coll Surg Engl* 65: 139-144
  10. De Santis D., Dinauer D., Duke J., Erlich H. A., Holcomb C. L., et al. (2013) 16(th) IHIW : review of HLA typing by NGS. *Int J Immunogenet* 40: 72-76
-

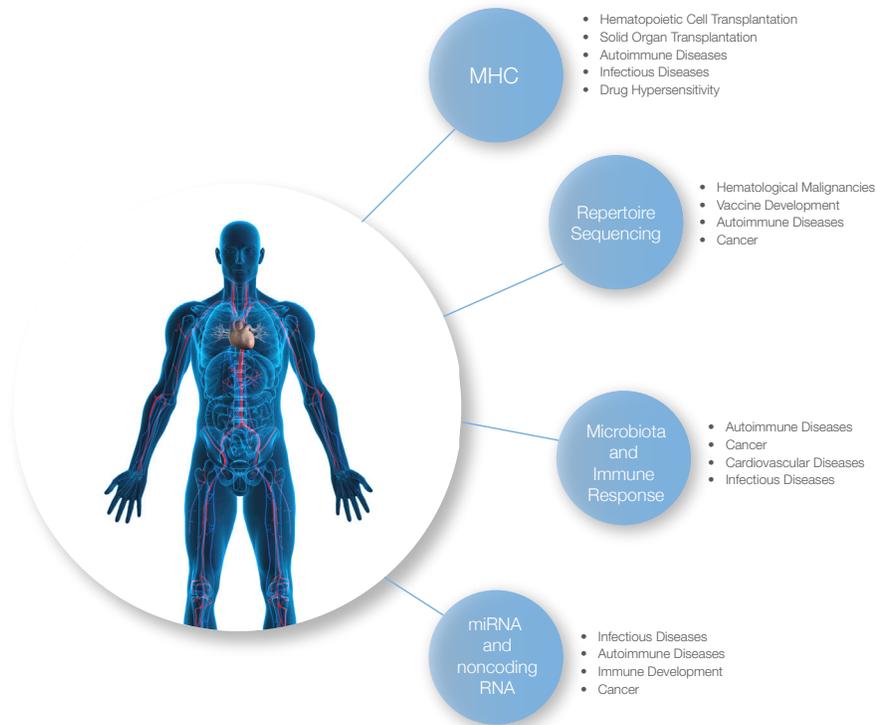

Next-generation sequencing has enabled researchers to make a significant impact in these applications (shown in blue). The list on the far right represents a subset of the human health and disease issues that can be addressed with these applications.

## Reviews

Georgiou G., Ippolito G. C., Beausang J., Busse C. E., Wardemann H., et al. (2014) The promise and challenge of high-throughput sequencing of the antibody repertoire. *Nat Biotechnol* 32: 158-168

Ansel K. M. (2013) RNA regulation of the immune system. *Immunol Rev* 253: 5-11

Bronevetsky Y. and Ansel K. M. (2013) Regulation of miRNA biogenesis and turnover in the immune system. *Immunol Rev* 253: 304-316

De Santis D., Dinauer D., Duke J., Erlich H. A., Holcomb C. L., et al. (2013) 16(th) IHIW : review of HLA typing by NGS. *Int J Immunogenet* 40: 72-76

Finn J. A. and Crowe J. E., Jr. (2013) Impact of new sequencing technologies on studies of the human B cell repertoire. *Curr Opin Immunol* 25: 613-618

Robins H. (2013) Immunosequencing: applications of immune repertoire deep sequencing. *Curr Opin Immunol* 25: 646-652

Shay T. and Kang J. (2013) Immunological Genome Project and systems immunology. *Trends Immunol* 34: 602-609

Warren E. H., Matsen F. A. t. and Chou J. (2013) High-throughput sequencing of B- and T-lymphocyte antigen receptors in hematology. *Blood* 122: 19-22

Woodsworth D. J., Castellarin M. and Holt R. A. (2013) Sequence analysis of T-cell repertoires in health and disease. *Genome Med* 5: 98


# ADAPTIVE IMMUNITY

## Repertoire Sequencing of Lymphocyte Receptors


The B and T-cell lymphocytes constitute the adaptive branch of the immune system, which is capable of identifying a vast range of antigens. This diverse repertoire of recognition elements is created through unique arrangements of immunoglobulin molecules on B-cell and T-cell receptors. Successful recognition of antigens triggers both an effector immune response as well as a memory response. An effector response includes CD8+ T-cells that eliminate cells with foreign antigens and CD4+ T-cells that differentiate into several different kinds of effector cells, including those that can further activate macrophages, cytotoxic T-cells, and B cells.<sup>11,12</sup> The B-cell effector response involves plasma cells that secrete antibodies capable of neutralizing or eliminating a foreign agent.<sup>13</sup> The memory response occurs when B and T-cells are activated by exposure to a foreign antigen.<sup>14</sup> Activation of these cells results in proliferation and preservation of the specific antigen receptor, such that secondary exposure to the foreign agent results in a robust immune response.<sup>15</sup>

In comparison to somatic cells, B and T-cell lymphocytes are unique in that their development and maturation are determined by DNA sequences that are not encoded in the germline. Instead, during the maturation process, these cells undergo rearrangement of the variable (V), diversity (D) and joining (J) gene segments in order to create a unique sequence that can encode an exclusive receptor structure in the heavy immunoglobulin chain of B cells, the  $\beta$  chain of  $\alpha\beta$  T-cell receptors, and the  $\delta$  chain of  $\gamma\delta$  T-cell receptors.

11. Litman G. W., Rast J. P. and Fugmann S. D. (2010) The origins of vertebrate adaptive immunity. *Nat Rev Immunol* 10: 543-553
12. Zhou L., Chong M. M. and Littman D. R. (2009) Plasticity of CD4+ T cell lineage differentiation. *Immunity* 30: 646-655
13. Lund F. E. and Randall T. D. (2010) Effector and regulatory B cells: modulators of CD4+ T cell immunity. *Nat Rev Immunol* 10: 236-247
14. Tokoyoda K., Hauser A. E., Nakayama T. and Radbruch A. (2010) Organization of immunological memory by bone marrow stroma. *Nat Rev Immunol* 10: 193-200
15. Ahmed R. and Gray D. (1996) Immunological memory and protective immunity: understanding their relation. *Science* 272: 54-60
16. Woodsworth D. J., Castellarin M. and Holt R. A. (2013) Sequence analysis of T-cell repertoires in health and disease. *Genome Med* 5: 98

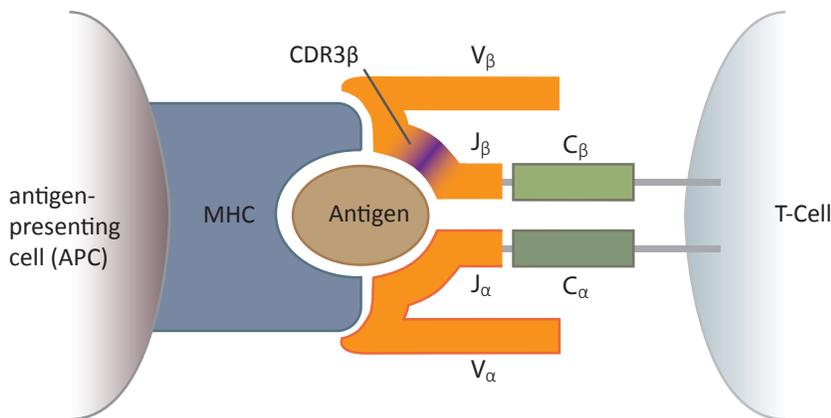


T-cell receptor-antigen-peptide-MHC interaction and T-cell receptor (TCR) gene recombination. (a) The antigen-presenting cell presents the peptide antigen bound to the major histocompatibility complex (MHC). The TCR (orange) binds to both the antigen and MHC. If the binding avidity is sufficiently high the T-cell is activated. The complementarity determining region 3 (CDR3) domain is shown in purple.<sup>16</sup>

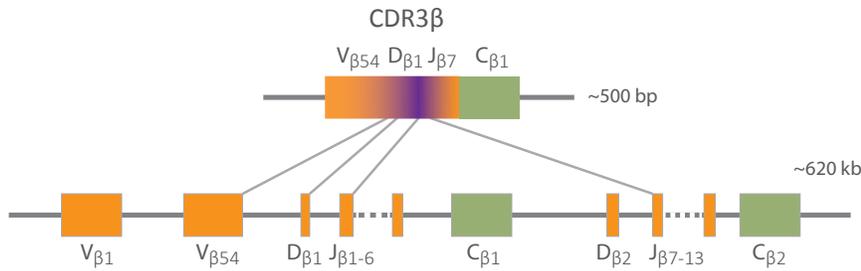

Simplified representation of TCR- $\beta$  VDJ gene recombination resulting in TCR diversity. The TCR- $\beta$  locus is located on chromosome 7 and is approximately 620 kb in length. Initially one of the two D regions is joined with one of 13 J regions (both randomly selected), followed by joining of the DJ region to one of more than 50 V regions (also randomly selected), yielding a final VDJ region that is approximately 500 bps in length. The mechanism by which gene segments are joined also introduces bp variability, which together with the combinatorial selection of these segments results in TCR diversity. A completely analogous process occurs for the TCR  $\alpha$  chain, without the D gene segment included.

17. Georgiou G., Ippolito G. C., Beausang J., Busse C. E., Wardemann H., et al. (2014) The promise and challenge of high-throughput sequencing of the antibody repertoire. *Nat Biotechnol* 32: 158-168

VDJ rearrangement in a B-cell generates the variable heavy chain of the immunoglobulin molecule. This immunoglobulin molecule is expressed on the surface of B cells and can also be freely secreted as an antibody.

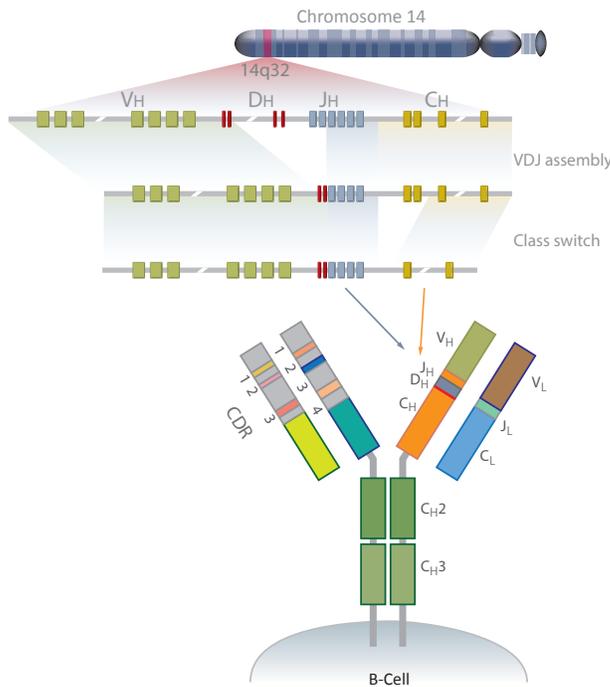

The primary antibody heavy chain repertoire is created predominantly by the somatic recombination of variable (V), diversity (D) and joining (J) gene segments. Nontemplated nucleotides (indicated in red) can also be added. The antigen-binding site of a heavy chain is formed by the juxtaposition of the hypervariable complementarity-determining regions (CDR-H1, H2 and H3) and the framework 3 region (FR3). After productive IgH rearrangement, recombination of the light chain (IgL) ensues, and the heterodimeric pairing of H and L chains forms the complete antibody of the IgM isotype that is expressed on the surface of a newly formed immature B cell.<sup>17</sup>

In addition to the aforementioned combinatorial diversity, made possible by the rearrangement of V, D, and J gene segments, splice variants contribute to this diversity via template-independent insertion and deletion of nucleotides at the V-D, D-J, and V-J splice junctions.<sup>18</sup> The size of the repertoire is further increased in B-cell receptors by somatic hypermutations (SHM) of B-cell receptor genes during affinity maturation after initial antigen encounter.

This combinatorial mechanism has the potential to generate more than  $10^{18}$  unique T-cell receptors in humans and a much more diverse B-cell repertoire.<sup>19,20</sup> The entire human VDJ region is has been estimated to range from 300 to 400 nucleotides in length, which makes read length a critical parameter in high-throughput sequencing.<sup>21,22,26</sup>

Complementarity determining regions (CDRs) are regions within antibodies or T cell receptors that complement an antigen's shape. Of the three complementarity-determining regions, CDR3 is the most variable locus and the most critical determinant of antigenic specificity.<sup>23-26</sup> The CDR3 region in the majority of rearranged functional TCR  $\beta$  and immunoglobulin heavy chains has a length ranging from 66<sup>21,27</sup> to 90<sup>22</sup> bps. Therefore, a sequencing depth of  $1 \times 10^9$  successfully characterizes the entire B and T-cell repertoire.<sup>28</sup>

Repertoire sequencing has applications in characterizing reconstitution of B and T-cell repertoires after hematopoietic stem cell transplantation, tracking lymphocytes in hematological malignancies, assessing vaccine efficacy, identifying lymphocyte repertoire variants associated with autoimmune diseases, and in identifying lymphocyte receptor variants in cancers such as colorectal cancer.<sup>29,30</sup>

18. Martinez N. M. and Lynch K. W. (2013) Control of alternative splicing in immune responses: many regulators, many predictions, much still to learn. *Immunol Rev* 253: 216-236
19. Vahedi G., Takahashi H., Nakayama S., Sun H. W., Sartorelli V., et al. (2012) STATs Shape the Active Enhancer Landscape of T Cell Populations. *Cell* 151: 981-993
20. Venturi V., Price D. A., Douek D. C. and Davenport M. P. (2008) The molecular basis for public T-cell responses? *Nat Rev Immunol* 8: 231-238
21. Rocha P. P., Micsinai M., Kim J. R., Hewitt S. L., Souza P. P., et al. (2012) Close proximity to Igh is a contributing factor to AID-mediated translocations. *Mol Cell* 47: 873-885
22. Benichou J., Ben-Hamo R., Louzoun Y. and Efroni S. (2012) Rep-Seq: uncovering the immunological repertoire through next-generation sequencing. *Immunology* 135: 183-191
23. Warren R. L., Freeman J. D., Zeng T., Choe G., Munro S., et al. (2011) Exhaustive T-cell repertoire sequencing of human peripheral blood samples reveals signatures of antigen selection and a directly measured repertoire size of at least 1 million clonotypes. *Genome Res* 21: 790-797
24. Robins H. S., Campregher P. V., Srivastava S. K., Wacher A., Turtle C. J., et al. (2009) Comprehensive assessment of T-cell receptor beta-chain diversity in alphabeta T cells. *Blood* 114: 4099-4107
25. Robins H. S., Srivastava S. K., Campregher P. V., Turtle C. J., Andriessen J., et al. (2010) Overlap and effective size of the human CD8+ T cell receptor repertoire. *Sci Transl Med* 2: 47ra64
26. Wang C., Sanders C. M., Yang Q., Schroeder H. W., Jr., Wang E., et al. (2010) High throughput sequencing reveals a complex pattern of dynamic interrelationships among human T cell subsets. *Proc Natl Acad Sci U S A* 107: 1518-1523
27. Larimore K., McCormick M. W., Robins H. S. and Greenberg P. D. (2012) Shaping of human germline IgH repertoires revealed by deep sequencing. *J Immunol* 189: 3221-3230
28. Warren E. H., Matsen F. A. t. and Chou J. (2013) High-throughput sequencing of B- and T-lymphocyte antigen receptors in hematology. *Blood* 122: 19-22
29. Benichou J., Ben-Hamo R., Louzoun Y. and Efroni S. (2012) Rep-Seq: uncovering the immunological repertoire through next-generation sequencing. *Immunology* 135: 183-191
30. Sherwood A. M., Emerson R. O., Scherer D., Habermann N., Buck K., et al. (2013) Tumor-infiltrating lymphocytes in colorectal tumors display a diversity of T cell receptor sequences that differ from the T cells in adjacent mucosal tissue. *Cancer Immunol Immunother* 62: 1453-1461

## Reviews

- Georgiou G., Ippolito G. C., Beausang J., Busse C. E., Wardemann H., et al. (2014) The promise and challenge of high-throughput sequencing of the antibody repertoire. *Nat Biotechnol* 32: 158-168
- Finn J. A. and Crowe J. E., Jr. (2013) Impact of new sequencing technologies on studies of the human B cell repertoire. *Curr Opin Immunol* 25: 613-618
- Robins H. (2013) Immunosequencing: applications of immune repertoire deep sequencing. *Curr Opin Immunol* 25: 646-652
- Warren E. H., Matsen F. A. t. and Chou J. (2013) High-throughput sequencing of B- and T-lymphocyte antigen receptors in hematology. *Blood* 122: 19-22
- Woodsworth D. J., Castellarin M. and Holt R. A. (2013) Sequence analysis of T-cell repertoires in health and disease. *Genome Med* 5: 98
- Benichou J., Ben-Hamo R., Louzoun Y. and Efroni S. (2012) Rep-Seq: uncovering the immunological repertoire through next-generation sequencing. *Immunology* 135: 183-191

## References

- Emerson R., Sherwood A., Desmarais C., Malhotra S., Phippard D., et al. (2013) Estimating the ratio of CD4+ to CD8+ T cells using high-throughput sequence data. *J Immunol Methods* 391: 14-21**  
The authors identify sequence features in the variable CDR3 region of the rearranged T cell receptor gene that distinguish CD4+ from CD8+ T cells. These features include variable gene usage and CDR3 region length. They estimate that as few as 1000 T cell receptor sequences are needed to accurately estimate the proportion of CD4+ and CD8+ T cells.  
Illumina Technology: HiSeq 2000
- Putintseva E. V., Britanova O. V., Staroverov D. B., Merzlyak E. M., Turchaninova M. A., et al. (2013) Mother and child T cell receptor repertoires: deep profiling study. *Front Immunol* 4: 463**  
The authors performed comparative analysis of these TCR repertoires of 3 mothers and 6 children. Thymic selection shapes the initial output of the TCR recombination machinery in both related and unrelated pairs, with minor effect from inherited differences. TCR profiling using characteristic TCR beta CDR3 variants as clonal identifiers also showed that mature T cells, transferred across the placenta during pregnancy, can expand and persist as functional microchimeric clones in their new host.  
Illumina Technology: HiSeq 2000
- Medvedovic J., Ebert A., Tagoh H., Tamir I. M., Schwickert T. A., et al. (2013) Flexible long-range loops in the VH gene region of the Igh locus facilitate the generation of a diverse antibody repertoire. *Immunity* 39: 229-244
- Meier J., Roberts C., Avent K., Hazlett A., Berrie J., et al. (2013) Fractal organization of the human T cell repertoire in health and after stem cell transplantation. *Biol Blood Marrow Transplant* 19: 366-377
- Genolet R., Stevenson B. J., Farinelli L., Osteras M. and Luescher I. F. (2012) Highly diverse TCRalpha chain repertoire of pre-immune CD8(+) T cells reveals new insights in gene recombination. *EMBO J* 31: 1666-1678
- Larimore K., McCormick M. W., Robins H. S. and Greenberg P. D. (2012) Shaping of human germline IgH repertoires revealed by deep sequencing. *J Immunol* 189: 3221-3230
- Wu D., Sherwood A., Fromm J. R., Winter S. S., Dunsmore K. P., et al. (2012) High-throughput sequencing detects minimal residual disease in acute T lymphoblastic leukemia. *Sci Transl Med* 4: 134ra163v

## Notes on Experimental Design

The primary challenges in CDR3 sequencing are the accumulation of PCR errors, sequencing errors, and ratio bias. These factors can result in the generation of false TCR diversity (artificial diversity) with the resultant inability to interpret sequence information accurately. For example, generating more sequencing reads may result in the expansion of erroneous sequence variants with one, two, or more mismatches. These amplified errors may be interpreted as evidence for sequence diversity.<sup>31,32</sup>

To address artificial diversity, previous studies suggest withdrawal of the low-abundance CDR3 variants that differ from the high-abundance variants by a single nucleotide mismatch or blind elimination of low-abundance sequence variants that comprise a total of 4% of all sequencing reads. It has been shown that this approach can result in up to 50% loss of sequencing reads and an even greater loss on non-Illumina platforms.<sup>33</sup>

To eliminate PCR and sequencing errors the following recommendations have been made in the literature.<sup>34</sup>

- From each sequencing read the CDR3 is extracted by aligning each sequence to the set of genomic VDJ segments from the IMGT/GENE-DB database. Low-quality nucleotides for VDJ segments are treated as allowable mismatches.
- Mapping low-quality reads. High-quality sequences at each nucleotide position within CDR3 form “core clonotypes.” These are merged with low-quality sequencing reads that have  $\leq 3$  low-quality nucleotides.
- Correcting PCR errors. Given that TCRs do not undergo somatic hypermutation, nucleotide mismatches with the V, D, or J segments of CDR3 can only arise from PCR and sequencing errors. Low-abundant core clonotypes are merged with the more abundant (at least 5-fold more abundant) core clonotypes that differ by no more than 3 nucleotides.

mRNA is the preferable starting material for TCR profiling.<sup>35</sup>

- T-cell contains multiple copies of RNA molecules that encode beta and alpha chains. These copies widen the bottleneck between the sampled T cells and the final TCR amplicon.
- Given that genomic DNA requires that the entire sample be amplified to compute TCR repertoire, this becomes technically challenging when studying sizable populations of T-cells, which would require unreasonably large aliquot volumes.

- 
31. Warren R. L., Freeman J. D., Zeng T., Choe G., Munro S., et al. (2011) Exhaustive T-cell repertoire sequencing of human peripheral blood samples reveals signatures of antigen selection and a directly measured repertoire size of at least 1 million clonotypes. *Genome Res* 21: 790-797
  32. Nguyen P., Ma J., Pei D., Obert C., Cheng C., et al. (2011) Identification of errors introduced during high throughput sequencing of the T cell receptor repertoire. *BMC Genomics* 12: 106
  33. Bolotin D. A., Mamedov I. Z., Britanova O. V., Zvyagin I. V., Shagin D., et al. (2012) Next generation sequencing for TCR repertoire profiling: platform-specific features and correction algorithms. *Eur J Immunol* 42: 3073-3083
  34. Bolotin D. A., Mamedov I. Z., Britanova O. V., Zvyagin I. V., Shagin D., et al. (2012) Next generation sequencing for TCR repertoire profiling: platform-specific features and correction algorithms. *Eur J Immunol* 42: 3073-3083
  35. Bolotin D. A., Mamedov I. Z., Britanova O. V., Zvyagin I. V., Shagin D., et al. (2012) Next generation sequencing for TCR repertoire profiling: platform-specific features and correction algorithms. *Eur J Immunol* 42: 3073-3083
-


## Single-Cell Repertoire Sequencing

Analysis of the immunoglobulin variable region and T-cell receptor repertoires is of fundamental importance for our understanding of adaptive immunity in health and disease.<sup>36</sup> However, the vast majority of repertoire studies yield data on only one of the two chains of immune receptors and thus cannot provide information about the identity of native receptor pairs encoded by single B or T-cell lymphocytes.<sup>37,38</sup>

Phage and yeast display technologies,<sup>39-41</sup> although efficient for isolation of antigen-specific antibodies, rely on random pairing and do not provide information on the native pairs of chains. Methods which involve growing cultures of lymphocyte clones,<sup>42</sup> or sorting of narrow antigen-specific populations of T cells<sup>43</sup> or B cells<sup>44</sup> are limited by the number of clones that can be identified, as well as by the complexity of biological samples.

New approaches to this problem take advantage of the sensitivity of next-generation sequencing to sequence single cells and identify multiple native TCR chain pairs in a single experiment.<sup>45</sup>

- 
36. Miles J. J., Douek D. C. and Price D. A. (2011) Bias in the alphabeta T-cell repertoire: implications for disease pathogenesis and vaccination. *Immunol Cell Biol* 89: 375-387
  37. Wu X., Zhou T., Zhu J., Zhang B., Georgiev I., et al. (2011) Focused evolution of HIV-1 neutralizing antibodies revealed by structures and deep sequencing. *Science* 333: 1593-1602
  38. Fischer N. (2011) Sequencing antibody repertoires: the next generation. *MAbs* 3: 17-20
  39. Bowley D. R., Jones T. M., Burton D. R. and Lerner R. A. (2009) Libraries against libraries for combinatorial selection of replicating antigen-antibody pairs. *Proc Natl Acad Sci U S A* 106: 1380-1385
  40. Marks J. D., Hoogenboom H. R., Bonnett T. P., McCafferty J., Griffiths A. D., et al. (1991) By-passing immunization. Human antibodies from V-gene libraries displayed on phage. *J Mol Biol* 222: 581-597
  41. Hoogenboom H. R., Griffiths A. D., Johnson K. S., Chiswell D. J., Hudson P., et al. (1991) Multi-subunit proteins on the surface of filamentous phage: methodologies for displaying antibody (Fab) heavy and light chains. *Nucleic Acids Res* 19: 4133-4137
  42. Lagerkvist A. C., Furebring C. and Borrebaeck C. A. (1995) Single, antigen-specific B cells used to generate Fab fragments using CD40-mediated amplification or direct PCR cloning. *Biotechniques* 18: 862-869
  43. Trautmann L., Rimbart M., Echasserieau K., Saulquin X., Neveu B., et al. (2005) Selection of T cell clones expressing high-affinity public TCRs within Human cytomegalovirus-specific CD8 T cell responses. *J Immunol* 175: 6123-6132
  44. Franz B., May K. F., Jr., Dranoff G. and Wucherpfennig K. (2011) Ex vivo characterization and isolation of rare memory B cells with antigen tetramers. *Blood* 118: 348-357
  45. Turchaninova M. A., Britanova O. V., Bolotin D. A., Shugay M., Putintseva E. V., et al. (2013) Pairing of T-cell receptor chains via emulsion PCR. *Eur J Immunol* 43: 2507-2515
-

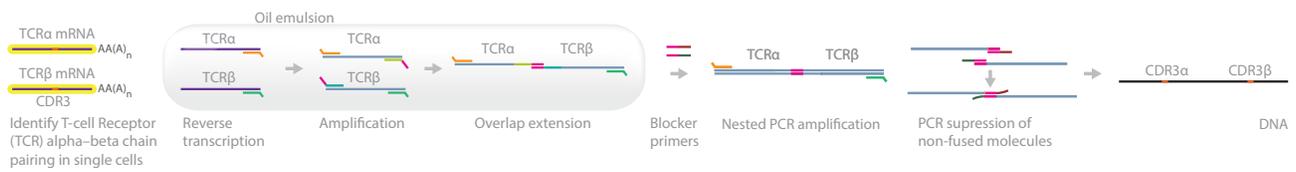

Cell-based emulsion RT-PCR technique for identifying TCR alpha-beta chain pairing. Released TCR alpha and beta mRNAs are reverse-transcribed, amplified, and overlap extended within each droplet. Products are extracted from the emulsion and fused molecules of interest are selectively amplified. Non-fused molecules are suppressed with blocking primers.<sup>46</sup>

## References

Dekosky B. J., Ippolito G. C., Deschner R. P., Lavinder J. J., Wi ne Y., et al. (2013) High-throughput sequencing of the paired human immunoglobulin heavy and light chain repertoire. *Nat Biotechnol* 31: 166-169

Previously VH:VL pairing in B-cell repertoire diversity was lost during lysis of B-cell populations. Here the authors employed a method of single-cell mRNA capture, reverse transcription and amplification by emulsion VH:VL linkage RT PCR of these pairings. The linked pairings were sequenced to identify unique antibody clonotypes in healthy peripheral blood IgG+ B-cells, peripheral antigen-specific plasmablasts isolated after tetanus toxoid immunization, and memory B-cell responses following influenza vaccination.

Illumina Technology: MiSeq 2 x 250 bp

Han A., Newell E. W., Glanville J., Fernandez-Becker N., Khosla C., et al. (2013) Dietary gluten triggers concomitant activation of CD4+ and CD8+ alphabeta T cells and gammadelta T cells in celiac disease. *Proc Natl Acad Sci U S A* 110: 13073-13078

Celiac disease is an intestinal autoimmune disease caused by dietary gluten and gluten-specific CD4+ T-cell responses. Gluten exposure also induces the appearance of activated, gut-homing CD8+ αβ and γδ T cells in peripheral blood. Single-cell T-cell receptor sequence analysis indicates that both of these cell populations have highly focused Tcell receptor repertoires. Such a focused repertoire usually indicates that the induction is driven by an antigen.

Illumina Technology: MiSeq paired-end sequencing

Turchaninova M. A., Britanova O. V., Bolotin D. A., Shugay M., Putintseva E. V., et al. (2013) Pairing of T-cell receptor chains via emulsion PCR. *Eur J Immunol* 43: 2507-2515

The authors propose a single cell-based method to identify native pairs of alpha-beta T cell receptor (TCR) CDR3 chains within emulsion droplets by employing reverse-transcription of alpha and beta chain mRNA, PCR amplification, and subsequent fusion via overlap-extension. This PCR suppression technique resolves the issue of random overlap-extension of gene pairs that may create a high level of noise after the emulsion stage. The authors propose that this methodology can be applied to the identification of native pairs of variable heavy-light antibody chains.

Illumina Technology: MiSeq 2 x 150 bp

- 
46. Turchaninova M. A., Britanova O. V., Bolotin D. A., Shugay M., Putintseva E. V., et al. (2013) Pairing of T-cell receptor chains via emulsion PCR. *Eur J Immunol* 43: 2507-2515
-

# LYMPHOCYTE DEVELOPMENT

## T Cell Development

Multipotent or lymphoid-biased precursors enter the T cell developmental pathway in response to signals from the thymic microenvironment.<sup>47</sup> Studies have shown that Notch, which has been classically associated with embryonic cell development, is an important trigger in T-cell lineage commitment. Notch signaling in the thymus causes hematopoietic precursors to commit to the T cell fate, mobilizes a T-cell gene expression program that prepares the cells for T-cell antigen receptor (TCR), TCR-based repertoire selection, and ultimately, prepares them for functional roles as immune effectors.<sup>48</sup>

Many questions remain regarding the molecular mechanisms of this commitment. For example, precursors entering the thymus display regulatory genes that are either expressed or inducible, however upon commitment these genes are not only repressed but also irreversibly silenced.<sup>49</sup>

Other questions relate to multiple regulatory requirements for successful deployment of the T-cell program. For example, there is a need to elucidate the functional role of additional transcription factors, including E2A, and HEB, TCF-1 and LEF-1, GATA-3, Myb, Runx1, Ikaros, and Gfi1.<sup>50</sup>


### Reviews

Martinez N. M. and Lynch K. W. (2013) Control of alternative splicing in immune responses: many regulators, many predictions, much still to learn. *Immunol Rev* 253: 216-236

Pagani M., Rossetti G., Panzeri I., de Candia P., Bonnal R. J., et al. (2013) Role of microRNAs and long-non-coding RNAs in CD4(+) T-cell differentiation. *Immunol Rev* 253: 82-96

Rothenberg E. V. (2012) Transcriptional drivers of the T-cell lineage program. *Curr Opin Immunol* 24: 132-138

- 
47. Petrie H. T. and Zuniga-Pflucker J. C. (2007) Zoned out: functional mapping of stromal signaling microenvironments in the thymus. *Annu Rev Immunol* 25: 649-679
  48. Zhang J. A., Mortazavi A., Williams B. A., Wold B. J. and Rothenberg E. V. (2012) Dynamic Transformations of Genome-wide Epigenetic Marking and Transcriptional Control Establish T Cell Identity. *Cell* 149: 467-482
  49. Vigano M. A., Ivanek R., Balwierz P., Berninger P., van Nimwegen E., et al. (2013) An epigenetic profile of early T-cell development from multipotent progenitors to committed T-cell descendants. *Eur J Immunol*
  50. Rothenberg E. V., Moore J. E. and Yui M. A. (2008) Launching the T-cell-lineage developmental programme. *Nat Rev Immunol* 8: 9-21
- 

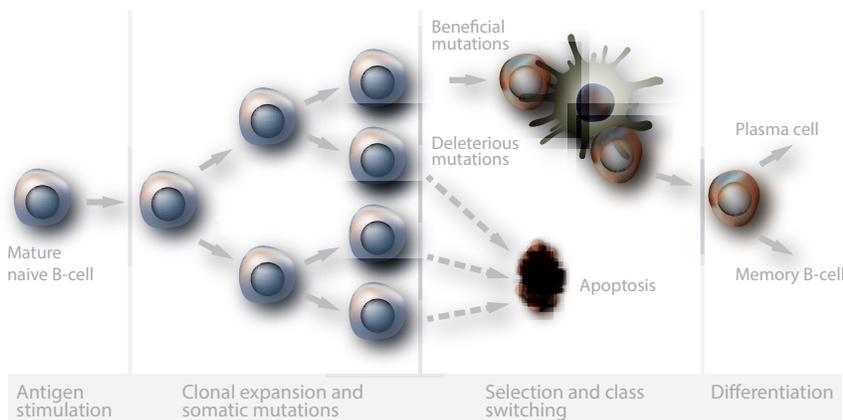

## References

Vigano M. A., Ivanek R., Balwierz P., Berninger P., van Nimwegen E., et al. (2013) An epigenetic profile of early T-cell development from multipotent progenitors to committed T-cell descendants. *Eur J Immunol*

In this study the authors analyzed and compared the gene expression profiles and the genome-wide histone modification marks H3K4me3 (H3 lysine 4 trimethylation) and H3K27me3 (H3 lysine 27 trimethylation) in T-cells cells at various stages of development. They observed global changes of gene expression and the epigenetic profile for H3K4me3 and H3K27me3 at promoters.

Illumina Technology: Genome Analyzer for ChIP-Seq analysis

Vahedi G., Takahashi H., Nakayama S., Sun H. W., Sartorelli V., et al. (2012) STATs Shape the Active Enhancer Landscape of T Cell Populations. *Cell* 151: 981-993

The authors mapped the activity of enhancer signatures using H3K4me1-high and p300-high regions in T helper 1 (Th1) and Th2 cells to interrogate active enhancer repertoires. They used RNA-Seq transcriptome profiling in T helper cells to identify the top 100 differentially expressed genes in each subset and did a comparative study between wild-type and STAT-deficient cells to assess STAT-dependency of positively regulated genes. They also used ChIP-Seq to demonstrate that STAT-deficient cells fail to fully recover the chromatin signature of STAT-dependent enhancers. Collectively, these findings suggest that STAT proteins play both a direct and indirect role in molding specialized enhancer architecture.

Illumina Technology: Genome Analyzer<sub>IIx</sub> for ChIP-Seq and HiSeq 2000 for RNA-Seq for 100 cycles (single read). RNA-Seq libraries were prepared using TruSeq sample prep kit

Zhang J. A., Mortazavi A., Williams B. A., Wold B. J. and Rothenberg E. V. (2012) Dynamic Transformations of Genome-wide Epigenetic Marking and Transcriptional Control Establish T Cell Identity. *Cell* 149: 467-482

Notch pathway signaling prompts hematopoietic precursors to become committed to the T cell fate and prepares cells for TCR expression and TCR-based repertoire selection. The authors employed RNA-Seq and ChIP-Seq to identify the dynamic transformations in transcription and epigenetic marking that occur across the genome through five stages of T cell differentiation that span lineage commitment (FLDN1, FLDN2a, FLDN2b, ThyDN3, and ThyDP cells). They reported that the major genome-wide transcriptomic changes leading to T lineage identity occur during transition to the DN2b or DN3; thus, during commitment and  $\beta$ -selection stages. The authors also used ChIP-Seq to enrich DNA associated with three H3 modifications: H3K(9, 14)Ac, H3K4me2, and H3K27me3. They employed these histone-marking patterns on potential cis-regulatory elements to *in vivo* tracking data of GATA-3 and PU.1 (transcription factors with complementary roles in early T cell development) and demonstrated the functional relevance of these transcription factors at potential sites for different developmental cell subsets.

Illumina Technology: Genome Analyzer for RNA-Seq and ChIP-Seq

Malinge S., Thiollier C., Chlon T. M., Doré L. C., Diebold L., et al. (2013) Ikaros inhibits megakaryopoiesis through functional interaction with GATA-1 and NOTCH signaling. *Blood* 121: 2440-2451

Genolet R., Stevenson B. J., Farinelli L., Osteras M. and Luescher I. F. (2012) Highly diverse TCRalpha chain repertoire of pre-immune CD8(+) T cells reveals new insights in gene recombination. *EMBO J* 31: 1666-1678


## B Cell Development

When mature B cells encounter an antigen they undergo a programmed DNA recombination event known as class switch recombination (CSR), which alters the effector function of the antibody molecule. During class switching, one constant region gene (typically C $\mu$ ) is replaced with another (either C $\gamma$ 3, C $\gamma$ 1, C $\gamma$ 2b, C $\gamma$ 2a, C $\epsilon$ , or C $\alpha$ ) via the introduction of double strand breaks (DSBs) and subsequent deletion of intervening sequences.<sup>51-53</sup>

In B lymphocytes, V(D)J recombination, class switch recombination (CSR) and somatic hypermutation (SHM) produce obligate single and double-strand DNA break intermediates that can become substrates for translocations.<sup>54,55</sup> These rearrangements could trigger cancer development.<sup>56</sup> This is supported by the observation that genetic ablation of the enzymes that create DNA lesions during V(D)J recombination (RAGs) or CSR and SHM (AID) has a significant protective effect on B-cell transformation.<sup>50,57</sup>

Nuclear architecture is another potential contributor to the incidence of chromosomal translocations.<sup>59</sup> Spatial organization of the genome is compartmentalized into chromosome territories as well as transcriptionally active and silent sub nuclear environments.<sup>60-63</sup> These compartments are believed to impact the frequency with which genes from different chromosomes interact and recombine.<sup>64</sup>

- 
51. Benichou J., Ben-Hamo R., Louzoun Y. and Efroni S. (2012) Rep-Seq: uncovering the immunological repertoire through next-generation sequencing. *Immunology* 135: 183-191
  52. Rocha P. P., Micsinai M., Kim J. R., Hewitt S. L., Souza P. P., et al. (2012) Close proximity to Igh is a contributing factor to AID-mediated translocations. *Mol Cell* 47: 873-885
  53. Larimore K., McCormick M. W., Robins H. S. and Greenberg P. D. (2012) Shaping of human germline IgH repertoires revealed by deep sequencing. *J Immunol* 189: 3221-3230
  54. Nussenzweig A. and Nussenzweig M. C. (2010) Origin of chromosomal translocations in lymphoid cancer. *Cell* 141: 27-38
  55. Tsai A. G., Lu H., Raghavan S. C., Muschen M., Hsieh C. L., et al. (2008) Human chromosomal translocations at CpG sites and a theoretical basis for their lineage and stage specificity. *Cell* 135: 1130-1142
  56. Mitelman F., Johansson B. and Mertens F. (2007) The impact of translocations and gene fusions on cancer causation. *Nat Rev Cancer* 7: 233-245
  57. Zhang Y., Gostissa M., Hildebrand D. G., Becker M. S., Boboila C., et al. (2010) The role of mechanistic factors in promoting chromosomal translocations found in lymphoid and other cancers. *Adv Immunol* 106: 93-133
  58. Klein I. A., Resch W., Jankovic M., Oliveira T., Yamane A., et al. (2011) Translocation-capture sequencing reveals the extent and nature of chromosomal rearrangements in B lymphocytes. *Cell* 147: 95-106
  59. Hakim O., Resch W., Yamane A., Klein I., Kieffer-Kwon K. R., et al. (2012) DNA damage defines sites of recurrent chromosomal translocations in B lymphocytes. *Nature* 484: 69-74
  60. Cremer T. and Cremer M. (2010) Chromosome territories. *Cold Spring Harb Perspect Biol* 2: a003889
  61. Lieberman-Aiden E., van Berkum N. L., Williams L., Imakaev M., Ragoczy T., et al. (2009) Comprehensive mapping of long-range interactions reveals folding principles of the human genome. *Science* 326: 289-293
  62. Hakim O., Sung M. H. and Hager G. L. (2010) 3D shortcuts to gene regulation. *Curr Opin Cell Biol* 22: 305-313
  63. Chakalova L. and Fraser P. (2010) Organization of transcription. *Cold Spring Harb Perspect Biol* 2: a000729
  64. Klein I. A., Resch W., Jankovic M., Oliveira T., Yamane A., et al. (2011) Translocation-capture sequencing reveals the extent and nature of chromosomal rearrangements in B lymphocytes. *Cell* 147: 95-106
-

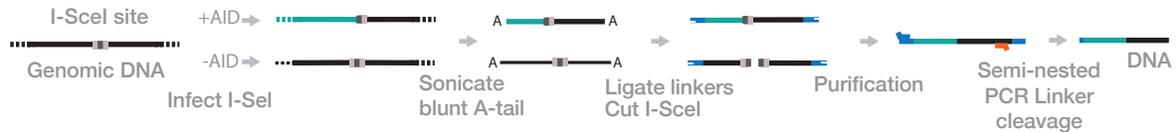

Translocation-capture sequencing (TC-Seq) is a method developed to study chromosomal rearrangements and translocations. In this method, cells are infected with retrovirus expressing I-SceI sites in cells with and without activation-induced cytidine deaminase (AICDA or AID) protein. Genomic DNA from cells is sonicated, linker-ligated, purified, and amplified via semi-nested LM-PCR. The linker is then cleaved and the DNA is sequenced. Any AID-dependent chromosomal rearrangement will be amplified by LM-PCR, while AID-independent translocations will be discarded.

## Reviews

Robbiani D. F. and Nussenzweig M. C. (2013) Chromosome translocation, B cell lymphoma, and activation-induced cytidine deaminase. *Annu Rev Pathol* 8: 79-103

## References

Holwerda S. J., van de Werken H. J., Ribeiro de Almeida C., Bergen I. M., de Bruijn M. J., et al. (2013) Allelic exclusion of the immunoglobulin heavy chain locus is independent of its nuclear localization in mature B cells. *Nucleic Acids Res* 41: 6905-6916

Chromatin conformation is one of many mechanisms for regulating gene expression. In developing B cells, the immunoglobulin heavy chain (IgH) locus undergoes a scheduled genomic rearrangement of the V, D, and J gene segments. In this study, an allele-specific chromosome conformation capture sequencing technique (4C-Seq) was applied to unambiguously follow the individual IgH alleles in mature B lymphocytes. The authors found that IgH adopts a lymphoid-specific nuclear location, and in mature B cells the distal VH regions of both IgH alleles position themselves away from active chromatin.

Illumina Technology: Genome Analyzer<sub>IIx</sub>, HiSeq 2000

Hakim O., Resch W., Yamane A., Klein I., Kieffer-Kwon K. R., et al. (2012) DNA damage defines sites of recurrent chromosomal translocations in B lymphocytes. *Nature* 484: 69-74

The authors performed chromosome conformation capture experiments followed by deep sequencing (4C-seq) to identify genomic regions that are in close spatial proximity to Igh and Myc, which are actively transcribed and targeted by AID. By comparing the chromosome conformation capture-on-chip (4C-seq) profiles to genome-wide epigenetic, transcription and TC-seq data sets, they concluded that Igh and Myc loci in AID-deficient peripheral B cells are more closely associated with epigenetically accessible genomic sites. They also employed a form of ChIP-Seq, termed RPA-seq, to measure the recruitment of replication protein A (RPA) in activated B cells. They used RPA recruitment as a proxy for AID-mediated DNA damage and demonstrated that the frequency of DNA damage directly accounts for the rate of translocation.

Illumina Technology: Genome Analyzer<sub>IIx</sub> for paired-end 4C-seq and ChIP-Seq

Rocha P. P., Micsinai M., Kim J. R., Hewitt S. L., Souza P. P., et al. (2012) Close proximity to Igh is a contributing factor to AID-mediated translocations. *Mol Cell* 47: 873-885

Approximately 95% of lymphomas are of B cell origin and many of these are attributed to aberrant rearrangements generated by activation-induced cytidine deaminase (AID) mediated breaks outside of the Igh locus. The authors used chromosome conformation capture-on-chip with massively parallel sequencing (4C-seq) to show that loci that have a significant interaction with Igh are specifically enriched for RNA Pol II, Spt5, H3K4me3, and AID. They employed a domain-centric approach for analyzing the 4C-seq data and found that chromosomal regions, which contact Igh at significant frequency, contain the vast majority (90%) of known sites identified as hotspot AID target genes. This study provides insight into the role of nuclear organization in AID targeting and the maintenance of genomic stability.

Illumina Technology: Genome Analyzer<sub>IIx</sub> for single-read 72-cycle run for 4C-seq and Illumina HiSeq 2000 to resequence a Cy1 library

Chaumeil J., Micsinai M., Ntziachristos P., Roth D. B., Aifantis I., et al. (2013) The RAG2 C-terminus and ATM protect genome integrity by controlling antigen receptor gene cleavage. *Nat Commun* 4: 2231

Jankovic M., Feldhahn N., Oliveira T. Y., Silva I. T., Kieffer-Kwon K. R., et al. (2013) 53BP1 alters the landscape of DNA rearrangements and suppresses AID-induced B cell lymphoma. *Mol Cell* 49: 623-631

Medvedovic J., Ebert A., Tagoh H., Tamir I. M., Schwickert T. A., et al. (2013) Flexible long-range loops in the VH gene region of the Igh locus facilitate the generation of a diverse antibody repertoire. *Immunity* 39: 229-244

Yamane A., Robbiani D. F., Resch W., Bothmer A., Nakahashi H., et al. (2013) RPA Accumulation during Class Switch Recombination Represents 5'-3' DNA-End Resection during the S-G2/M Phase of the Cell Cycle. *Cell Rep* 3: 138-147

Hakim O., Resch W., Yamane A., Klein I., Kieffer-Kwon K. R., et al. (2012) DNA damage defines sites of recurrent chromosomal translocations in B lymphocytes. *Nature* 484: 69-74

Kovalchuk A. L., Ansarah-Sobrinho C., Hakim O., Resch W., Tolarova H., et al. (2012) Mouse model of endemic Burkitt translocations reveals the long-range boundaries of Ig-mediated oncogene deregulation. *Proc Natl Acad Sci U S A* 109: 10972-10977


Rocha P. P., Micsinai M., Kim J. R., Hewitt S. L., Souza P. P., et al. (2012) Close proximity to Igh is a contributing factor to AID-mediated translocations. *Mol Cell* 47: 873-885

65. Noordermeer D., Branco M. R., Splinter E., Klous P., van Ijcken W., et al. (2008) Transcription and chromatin organization of a house-keeping gene cluster containing an integrated beta-globin locus control region. *PLoS Genet* 4: e1000016
66. Montavon T., Soshnikova N., Mascroz B., Joye E., Thevenet L., et al. (2011) A regulatory archipelago controls Hox genes transcription in digits. *Cell* 147: 1132-11453
67. Zhao Z., Tavosodana G., Sjolinder M., Gondor A., Mariano P., et al. (2006) Circular chromosome conformation capture (4C) uncovers extensive networks of epigenetically regulated intra- and interchromosomal interactions. *Nat Genet* 38: 1341-1347
68. Sajan S. A. and Hawkins R. D. (2012) Methods for identifying higher-order chromatin structure. *Annu Rev Genomics Hum Genet* 13: 59-826

## Notes on experimental design

4C-seq is the preferred chromosomal conformation capture technique when analyzing the DNA contact profile of individual genomic sites. This assay has been particularly useful for investigating the associations of specific genes with long-range regulatory elements.<sup>65,66</sup> (See **Chromatin Structure and Rearrangement** for more details).

4C is currently limited to the assessment of long-range contacts with larger regions elsewhere on the chromosome (in cis) or on other chromosomes (in trans). For example, local interactions (<50 kb distance) between a gene and its enhancer are not readily detected. Most 4C strategies use restriction enzymes with a 6-nucleotide recognition sequence, which cut once every few kilobases. This creates fragments that are much larger than the average regulatory sequences, which are no larger than several hundred bps. Increased resolution may depend on the use of more selective restriction enzymes that can generate shorter fragments, which can enable detection of de novo local regulatory interactions.

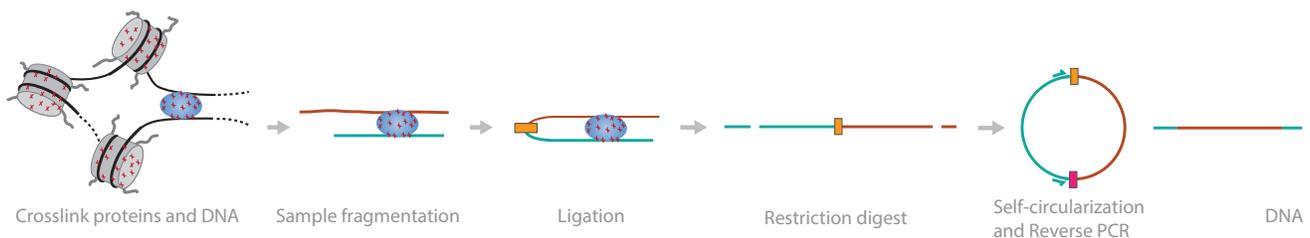

Circular chromatin conformation capture (4-C)<sup>67</sup>, allows the unbiased detection of all genomic regions that interact with a particular region of interest.<sup>68</sup> In this method, DNA-protein complexes are crosslinked using formaldehyde. The sample is fragmented, and the DNA is ligated and digested. The resulting DNA fragments self-circularize, followed by reverse PCR and sequencing. Deep sequencing provides base-pair resolution of ligated fragments.

# INNATE IMMUNITY

Innate immunity is the frontline of host defense. It prompts the rapid and local response against pathogens and is also important for the symbiotic partnership between the host and its microbiota. The innate immune system and adaptive immunity comprise the binary-classification of the immune response. Historically the distinction between these two branches of immunity relied on the consensus that innate immunity is nonspecific and lacks memory whereas the adaptive immunity is characterized by specific antigen recognition and subsequent memory response. However, new evidence that demonstrates innate immune features of B cells and T cells and adaptive immune properties of natural killer (NK) cells are now blurring this conventional binary-classification.<sup>69</sup>

Innate immunity involves the coordinated action of families of receptors, known as pattern-recognition receptors (PRRs) or microbial sensors that respond to a wide range of microorganisms through the detection of specific conserved microbial patterns or molecules.<sup>70-73</sup> This innate immune response is activated by specialized sets of receptors found on macrophages, mast cells, dendritic cells, natural killer cells, and polymorphonuclear leukocytes. These receptors include the membrane-bound Toll-like receptors (TLRs) and C-type lectin receptors (CLRs), and the cytosolic RIG-I-like receptors (RLRs), NOD-like receptors (NLRs) and other DNA sensors.<sup>74-78</sup> Furthermore, complement recognition molecules are circulating proteins that constitute the humoral (free in serum and body fluids) arm of innate immunity.<sup>79</sup> Following ligand binding, receptors induce the activation of distinct signaling pathways that involve effector molecules, such as interferons (IFNs) or antimicrobial peptides (AMPs), which are required for the eradication of pathogens or danger signals.

## Reviews

Quintana-Murci L. and Clark A. G. (2013) Population genetic tools for dissecting innate immunity in humans. *Nat Rev Immunol* 13: 280-293

- 
69. Lanier L. L. (2013) Shades of grey--the blurring view of innate and adaptive immunity. *Nat Rev Immunol* 13: 73-74
  70. Takeuchi O. and Akira S. (2010) Pattern recognition receptors and inflammation. *Cell* 140: 805-820
  71. Bae J., Ricciardi C. J., Esposito D., Komarnytsky S., Hu P., et al. (2014) Activation of Pattern Recognition Receptors in Brown Adipocytes Induce Inflammation and Suppress Uncoupling Protein 1 Expression and Mitochondrial Respiration. *Am J Physiol Cell Physiol*
  72. Walsh D., McCarthy J., O'Driscoll C. and Melgar S. (2013) Pattern recognition receptors--molecular orchestrators of inflammation in inflammatory bowel disease. *Cytokine Growth Factor Rev* 24: 91-104
  73. Fukata M. and Arditi M. (2013) The role of pattern recognition receptors in intestinal inflammation. *Mucosal Immunol* 6: 451-463
  74. Kawai T. and Akira S. (2011) Toll-like receptors and their crosstalk with other innate receptors in infection and immunity. *Immunity* 34: 637-650
  75. Maekawa T., Kufer T. A. and Schulze-Lefert P. (2011) NLR functions in plant and animal immune systems: so far and yet so close. *Nat Immunol* 12: 817-826
  76. Loo Y. M. and Gale M., Jr. (2011) Immune signaling by RIG-I-like receptors. *Immunity* 34: 680-692
  77. Schroder K. and Tschopp J. (2010) The inflammasomes. *Cell* 140: 821-832
  78. Sancho D. and Reis e Sousa C. (2012) Signaling by myeloid C-type lectin receptors in immunity and homeostasis. *Annu Rev Immunol* 30: 491-529
  79. Bottazzi B., Doni A., Garlanda C. and Mantovani A. (2010) An integrated view of humoral innate immunity: pentraxins as a paradigm. *Annu Rev Immunol* 28: 157-183
-

## References

Shalek A. K., Satija R., Adiconis X., Gertner R. S., Gaublot J. T., et al. (2013) Single-cell transcriptomics reveals bimodality in expression and splicing in immune cells. *Nature* 498: 236-240

The authors employed single-cell RNA sequencing to investigate heterogeneity in the response of mouse bone-marrow-derived dendritic cells (BMDCs) to lipopolysaccharides. Maturation of BMDCs can occur as a response to pathogen-derived ligands, such as lipopolysaccharides, which results in co-expression of defense cytokines. Here the authors identified a cluster of 137 highly variable, yet co-regulated, antiviral response genes. 100 out of the clusters' genes were bimodally expressed across single cells despite being strongly expressed at the population level. The authors stimulated and profiled BMDCs from interferon receptor knockout (*Irfn* <sup>-/-</sup>) mice and discovered markedly reduced expression in these cells for both *Stat2* and *Irf7*, as well other measured genes in the cluster. They concluded that interferon signaling is required for the induction of *Stat2* and *Irf7*, which subsequently act to induce variable antiviral cluster genes.

illumina Technology: HiSeq 2000 for single-cell RNA-Seq with average depth of 27 million read pairs

Cho P., Gelinas L., Corbett N. P., Tebbutt S. J., Turvey S. E., et al. (2013) Association of common single-nucleotide polymorphisms in innate immune genes with differences in TLR-induced cytokine production in neonates. *Genes Immun* 14: 199-211

Kemp C., Mueller S., Goto A., Barbier V., Paro S., et al. (2013) Broad RNA interference-mediated antiviral immunity and virus-specific inducible responses in *Drosophila*. *J Immunol* 190: 650-658

Lutay N., Ambite I., Gronberg Hernandez J., Rydstrom G., Ragnarsdottir B., et al. (2013) Bacterial control of host gene expression through RNA polymerase II. *J Clin Invest* 123: 2366-2379

Mikacenic C., Reiner A. P., Holden T. D., Nickerson D. A. and Wurfel M. M. (2013) Variation in the TLR10/TLR1/TLR6 locus is the major genetic determinant of interindividual difference in TLR1/2-mediated responses. *Genes Immun* 14: 52-57

Seed K. D., Lazinski D. W., Calderwood S. B. and Camilli A. (2013) A bacteriophage encodes its own CRISPR/Cas adaptive response to evade host innate immunity. *Nature* 494: 489-491

Wei Y., Chen R., Dimicoli S., Bueso-Ramos C., Neuberger D., et al. (2013) Global H3K4me3 genome mapping reveals alterations of innate immunity signaling and overexpression of JMJD3 in human myelodysplastic syndrome CD34+ cells. *Leukemia* 27: 2177-2186

Wie S. H., Du P., Luong T. Q., Rought S. E., Beliakova-Bethell N., et al. (2013) HIV Downregulates Interferon-Stimulated Genes in Primary Macrophages. *J Interferon Cytokine Res* 33: 90-95

Tsoi L. C., Spain S. L., Knight J., Ellinghaus E., Stuart P. E., et al. (2012) Identification of 15 new psoriasis susceptibility loci highlights the role of innate immunity. *Nat Genet* 44: 1341-1348

Yang G., Yang L., Zhao Z., Wang J. and Zhang X. (2012) Signature miRNAs Involved in the Innate Immunity of Invertebrates. *PLoS ONE* 7: e39015

# CANCER AND THE IMMUNE RESPONSE

The development from a normal hematopoietic cell to a cancerous cell involves a multistep process of clonal evolution driven by a series of somatic mutations. These mutations progressively transform the cell from normal growth to a precancerous state and finally a cancerous state, where all checkpoints designed to regulate cell growth have been surmounted.

Induction of malignant transformations appears to involve at least two distinct phases: initiation and promotion. Initiation involves changes in the genome but does not, in itself, lead to malignant transformation. Malignant transformation requires a secondary step, termed promotion. Promotion can occur during the aggressive cell division that follows the initiation phase, and results from the accumulation of new DNA alterations, typically affecting proto-oncogenes, tumor-suppressor genes or apoptotic genes, that result in unregulated cellular growth.

The ability of next-generation sequencing to detect mutations in rare clonal types, or cells, through deep sequencing makes it possible to study the role of immune effector functions in the pathogenesis of hematological malignancies. A notable example has been the influx of reports that implicate autoreactive T-cell clones in the pathogenesis of clonal stem cell disorders such as myelodysplastic syndromes (MDS) and aplastic anemia (AA).<sup>80</sup> These studies have been supported by the widely consolidated understanding that impairment of anti-tumor immunity, which is physiologically mediated by T-cells, can predispose the development of hematological malignancies. Collectively these T-cell repertoire studies and new reports that implicate immunoglobulin heavy chain rearrangements in clonal evolution of acute lymphoblastic leukemia have quickly become one of the most exciting research areas in hematology.<sup>81-83</sup>

## Reviews

Weaver W. M., Tseng P., Kunze A., Masaeli M., Chung A. J., et al. (2014) Advances in high-throughput single-cell microtechnologies. *Curr Opin Biotechnol* 25C: 114-123

Fozza C. and Longinotti M. (2013) T-cell receptor repertoire usage in hematologic malignancies. *Crit Rev Oncol Hematol* 86: 201-211

Warren E. H., Matsen F. A. t. and Chou J. (2013) High-throughput sequencing of B- and T-lymphocyte antigen receptors in hematology. *Blood* 122: 19-22

Gottgens B. (2012) Genome-scale technology driven advances to research into normal and malignant haematopoiesis. *Scientifica (Cairo)* 2012: 437956

For more information, see: [Cancer and the Immune System](#).


- 
80. Fozza, C., and Longinotti, M. (2013) T-cell receptor repertoire usage in hematologic malignancies. *Critical reviews in oncology/hematology* 86: 201–211
  81. Faham, M., Zheng, J., Moorhead, M., Carlton, V. E. H., Stow, P., et al. (2012) Deep-sequencing approach for minimal residual disease detection in acute lymphoblastic leukemia. *Blood* 120: 5173–5180
  82. Gawad, C., Pepin, F., Carlton, V. E. H., Klinger, M., Logan, A. C., et al. (2012) Massive evolution of the immunoglobulin heavy chain locus in children with B precursor acute lymphoblastic leukemia. *Blood* 120: 4407–4417
  83. Jan, M., Snyder, T. M., Corces-Zimmerman, M. R., Vyas, P., Weissman, I. L., et al. (2012) Clonal Evolution of Preleukemic Hematopoietic Stem Cells Precedes Human Acute Myeloid Leukemia. *Science Translational Medicine* 4: 149ra118–149ra118
-

# MICROBIOTA AND THE IMMUNE SYSTEM


Microbiota refers to the extraordinarily large and diverse reservoir of microorganisms that has a co-evolved relationship with the mammalian immune system. These complex microbial communities inhabit the body surfaces of virtually all vertebrates. The immune system plays an essential role in maintaining homeostasis with resident microbial communities to ensure that the mutualistic nature of the host-microbial relationships is sustained. The co-evolution of the vertebrate immune system has therefore been driven by the need to protect the host from pathogens and to foster complex microbial communities for their protective and metabolic benefits. Given that alterations in host-microbiota homeostasis have been implicated in viral infections,<sup>84,85</sup> autoimmune diseases,<sup>86</sup> cancer, metabolic diseases, and cardiovascular diseases, this is an exciting opportunity for researchers to examine the interactions between the microbiota and the host-immune response.

Next generation sequencing technologies have enabled researchers to define the construction of these microbiota by operationally defining polymorphisms of bacterial genes; especially those encoding the 16S ribosomal RNA sequences. Sequencing the human microbiome is now enabling researchers to examine the interactions between microbial communities and host immunity. This has illuminated the significant role that the immune system plays in mediating this homeostatic relationship.

- 
- 84. Sandler N. G., Koh C., Roque A., Eccleston J. L., Siegel R. B., et al. (2011) Host response to translocated microbial products predicts outcomes of patients with HBV or HCV infection. *Gastroenterology* 141: 1220-1230, 1230 e1221-1223
  - 85. Leung J. M., Davenport M., Wolff M. J., Wiens K. E., Abidi W. M., et al. (2014) IL-22-producing CD4+ cells are depleted in actively inflamed colitis tissue. *Mucosal Immunol* 7: 124-133
  - 86. McGuckin M. A., Eri R., Simms L. A., Florin T. H. and Radford-Smith G. (2009) Intestinal barrier dysfunction in inflammatory bowel diseases. *Inflamm Bowel Dis* 15: 100-113
- 

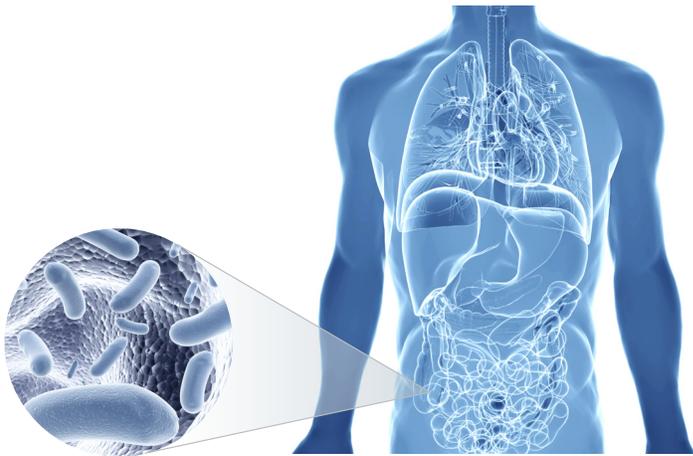


The human intestine harbors over 100 trillion microbes, which represent approximately 500 different species of bacteria.

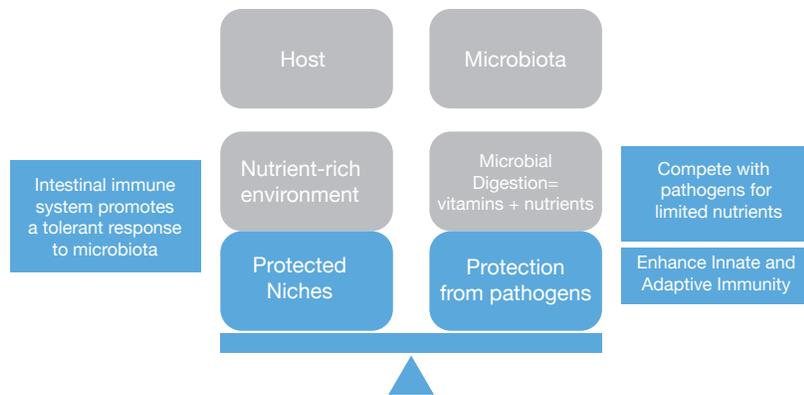

87. Kamada N., Seo S. U., Chen G. Y. and Nunez G. (2013) Role of the gut microbiota in immunity and inflammatory disease. *Nat Rev Immunol* 13: 321-335

The host and microbiota have co-evolved mutually beneficial outcomes and the immune system plays a critical role in preserving homeostasis. The host provides a nutrient-rich environment and protected niches for the microbiota. The microbiota provides the host with vitamins and nutrients as by-products of microbial digestion and protects the host from pathogens. The microbiota enhances the innate and adaptive immune response. Conversely, there is a need for the host to promote a tolerant immune response that enables the microbiota to inhabit the niches of the gut.

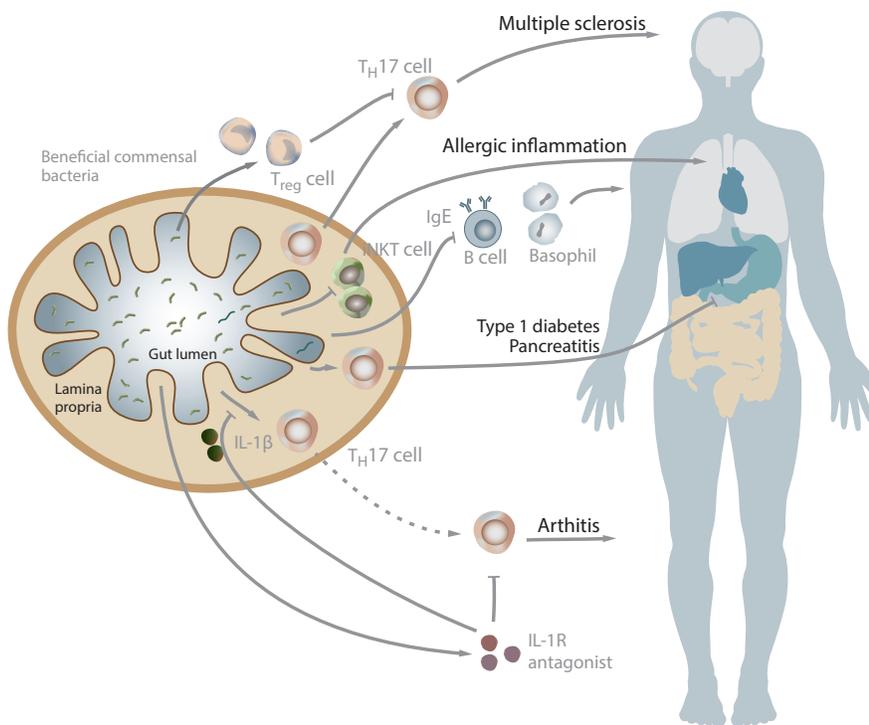

This network illustrates the various roles of the gut microbiota in extra-intestinal autoimmune diseases.<sup>87</sup>

## Reviews

Brown E. M., Sadarangani M. and Finlay B. B. (2013) The role of the immune system in governing host-microbe interactions in the intestine. *Nat Immunol* 14: 660-667

Castelino M., Eyre S., Upton M., Ho P. and Barton A. (2013) The bacterial skin microbiome in psoriatic arthritis, an unexplored link in pathogenesis: challenges and opportunities offered by recent technological advances. *Rheumatology (Oxford)*

Kamada N., Chen G. Y., Inohara N. and Nunez G. (2013) Control of pathogens and pathobionts by the gut microbiota. *Nat Immunol* 14: 685-690

Kamada N., Seo S. U., Chen G. Y. and Nunez G. (2013) Role of the gut microbiota in immunity and inflammatory disease. *Nat Rev Immunol* 13: 321-335

Mavrommatis B., Young G. R. and Kassiatis G. (2013) Counterpoise between the microbiome, host immune activation and pathology. *Curr Opin Immunol* 25: 456-462

Hooper L. V., Littman D. R. and Macpherson A. J. (2012) Interactions between the microbiota and the immune system. *Science* 336: 1268-1273

## References

[Leung J. M., Davenport M., Wolff M. J., Wiens K. E., Abidi W. M., et al. \(2014\) IL-22-producing CD4+ cells are depleted in actively inflamed colitis tissue. \*Mucosal Immunol\* 7: 124-133](#)

Changes to bacterial taxa may be associated with disease pathogenesis in inflammatory bowel disease. The authors performed deep sequencing analysis on the variable region 4 (V4) of bacterial 16S rRNA to investigate the mucosal microbiota communities in pinch biopsies from a cohort of ulcerative colitis (UC) patients. This deep sequencing study when combined with data generated by flow cytometry demonstrates that the depletion of Th22 cells, a subset of CD4+ helper T cells that produce IL-22, during active inflammation in UC patients is associated with reduced populations of Clostridiales and increased population of Proteobacteria among other specific alteration of the mucosal microbiota

Illumina Technology: MiSeq of 16S rRNA

[Markle J. G., Frank D. N., Mortin-Toth S., Robertson C. E., Feazel L. M., et al. \(2013\) Sex differences in the gut microbiome drive hormone-dependent regulation of autoimmunity. \*Science\* 339: 1084-1088](#)

Microbial factors, in particular the gut microbiota, are thought to influence susceptibility to type 1 diabetes. In the non-obese diabetic (NOD) mouse model of type 1 diabetes, female mice are significantly more susceptible to disease than males. This difference is not apparent under germ-free conditions, which indicates a role for microbes in type 1 diabetes susceptibility. The authors sequenced bacterial 16S rRNA libraries prepared from cecal contents to identify microbiome composition differences between sexes, through maturation and post-cecal transplantation. Transfer of cecal contents from male NOD mice to female NOD mice prior to disease onset protected against pancreatic islet inflammation, reduced autoantibody production, prevented the onset of diabetes, and was associated with increased testosterone in female mice. The protection conferred by M→F microbiome transfer was attenuated when androgen receptors were blocked. This study suggests that the microbiota may have a regulatory role on sex hormones and may influence autoimmune disease fate in individuals with high genetic risk.

Illumina Technology: MiSeq paired-end sequencing 16S rRNA

[Wang X., Lin Z., Gao L., Wang A., Wan Z., et al. \(2013\) Exome sequencing reveals a signal transducer and activator of transcription 1 \(STAT1\) mutation in a child with recalcitrant cutaneous fusariosis. \*J Allergy Clin Immunol\* 131: 1242-1243](#)

Fusarium species can cause papulopustular lesions with abscesses, ulcerations, and even necrosis in immunocompromised and immunocompetent patients. The authors report a case study of a 7 year old girl with treatment-resistant cutaneous fusariosis. They employed exome sequencing and identified a single heterozygous missense mutation in the signal transducer and activator of transcription 1 (STAT1) gene, which is the most likely genetic defect underlying the fusariosis in this patient.

Illumina Technology: HiSeq 2000 exome sequencing

Schwab C., Berry D., Rauch I., Rennisch I., Ramesmayer J., et al. (2014) Longitudinal study of murine microbiota activity and interactions with the host during acute inflammation and recovery. *ISME J*

Arthur J. C., Perez-Chanona E., Mühlbauer M., Tomkovich S., Uronis J. M., et al. (2012) Intestinal inflammation targets cancer-inducing activity of the microbiota. *Science* 338: 120-123

Caporaso J. G., Lauber C. L., Walters W. A., Berg-Lyons D., Huntley J., et al. (2012) Ultra-high-throughput microbial community analysis on the Illumina HiSeq and MiSeq platforms. *ISME J* 6: 1621-1624


Maslanik T., Tannura K., Mahaffey L., Loughridge A. B., Benninson L., et al. (2012) Commensal bacteria and MAMPs are necessary for stress-induced increases in IL-1 $\beta$  and IL-18 but not IL-6, IL-10 or MCP-1. *PLoS ONE* 7: e50636

# MAJOR HISTOCOMPATIBILITY COMPLEX

Although both T and B cells use surface receptors to recognize antigens, they accomplish this in two different ways. In contrast to antibodies or B-cell receptors, which can directly recognize antigens, T-cell receptors only recognize antigens that are presented on the surface of antigen presenting cells, such as dendritic cells and macrophages. These antigen peptides reside within the groove of a cell surface protein called the major histocompatibility complex (MHC) molecule.

In humans, the MHC locus is referred to as the Human Leukocyte Antigen (HLA) and encodes a collection of genes that span a contiguous 4 Mb region on the short arm of chromosome 6.<sup>88</sup> Moreover, the extended MHC, termed (xMHC), spans an even larger 7.6 Mb region comprising more than 400 annotated genes and pseudogenes.<sup>89</sup>

- 
- 88. de Bakker P. I. and Raychaudhuri S. (2012) Interrogating the major histocompatibility complex with high-throughput genomics. *Hum Mol Genet* 21: R29-36
  - 89. Horton R., Wilming L., Rand V., Lovering R. C., Bruford E. A., et al. (2004) Gene map of the extended human MHC. *Nat Rev Genet* 5: 889-899
- 

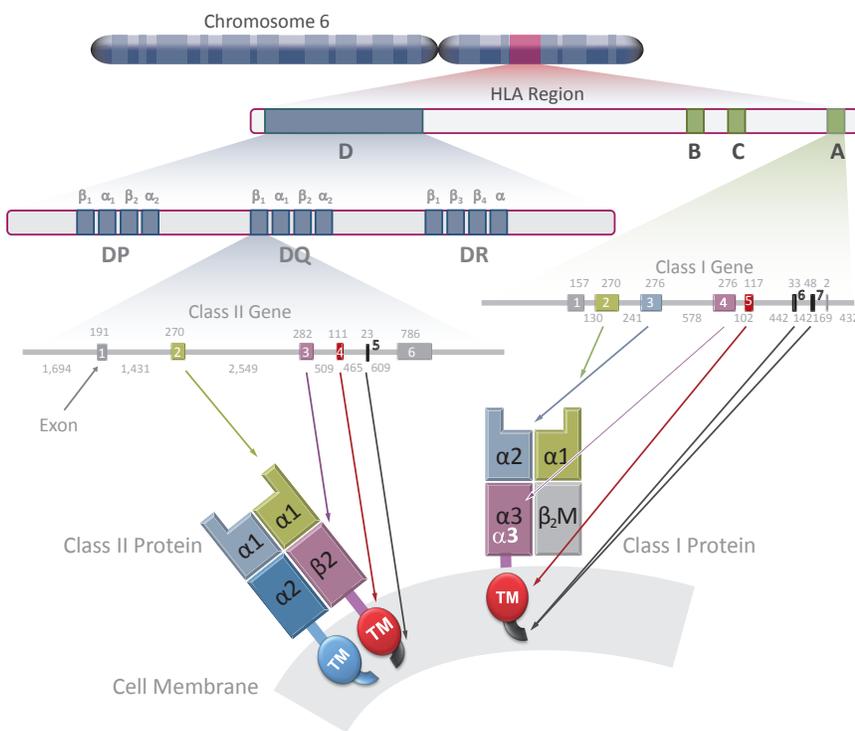

This is a map of the human MHC loci. The MHC class I genes are colored red, MHC class II genes are colored blue, and genes in MHC III are colored green. The 6 loci that are outlined encode the peptide binding sites of class I and II MHC molecules. These loci are routinely assessed in matching donors and recipients in hematopoietic cell and solid organ transplantation.

Six of the HLA genes (HLA-A, -B, -C, -DQA1, -DQB1 and -DRB1) are extremely polymorphic and constitute a set of important markers that are routinely employed in matching patients and donors for solid organ transplantation<sup>90</sup> and hematopoietic stem cell transplantation.<sup>91,92</sup> HLA genes also play an important role in infectious diseases (HIV, Hep C and CMV), autoimmune diseases (diabetes, rheumatoid arthritis, and celiac disease), and drug hypersensitivity.<sup>93</sup>

With conventional technologies, only the most polymorphic regions of HLA class I (exons 2 and 3) and II (exon 2), which encode the peptide binding sites, are assessed in clinical settings. Next-generation sequencing provides clinical researchers with the capability to sequence the entire gene, resulting in phase-resolved, unambiguous HLA typing. Clinical studies have shown that matching these regions of the 6 major HLA loci provides the best clinical outcomes with decreased incidence of rejection and graft versus host disease (GVHD) in solid organ and hematopoietic stem cell transplantation. However, even when these regions are matched, approximately 30% of recipients experience adverse events within 5 years.<sup>94</sup>

The source of these imperfect matches is unknown, but there are several possibilities. Adverse events may reflect mismatches in regions that lie outside of the regions that are currently analyzed. Given the high degree of polymorphisms, ambiguous combination of alleles may arise during HLA typing. These may result from cis/trans ambiguities or due to a particular allele combination being identical over the regions commonly analyzed. For example, in conventional sequencing both heterozygous alleles are coamplified and sequenced. Combination ambiguity can occur when two or more alleles share identical sequences in the targeted exons but exhibit differences in non-sequenced exons.<sup>95</sup> Next-generation sequencing provides phase information, which may significantly improve the analysis of the HLA cohort (See **Phase-Defined HLA Sequencing** for more details).

## Reviews

de Bakker P. I. and Raychaudhuri S. (2012) Interrogating the major histocompatibility complex with high-throughput genomics. *Hum Mol*

Erich H. (2012) HLA DNA typing: past, present, and future. *Tissue Antigens* 80: 1-11

- 
90. Eapen M., Rubinstein P., Zhang M. J., Stevens C., Kurtzberg J., et al. (2007) Outcomes of transplantation of unrelated donor umbilical cord blood and bone marrow in children with acute leukaemia: a comparison study. *Lancet* 369: 1947-1954
  91. Marks C. (1983) Immunobiological determinants in organ transplantation. *Ann R Coll Surg Engl* 65: 139-144
  92. Davies J. L., Kawaguchi Y., Bennett S. T., Copeman J. B., Cordell H. J., et al. (1994) A genome-wide search for human type 1 diabetes susceptibility genes. *Nature* 371: 130-136
  93. Mallal S., Nolan D., Witt C., Masel G., Martin A. M., et al. (2002) Association between presence of HLA-B\*5701, HLA-DR7, and HLA-DQ3 and hypersensitivity to HIV-1 reverse-transcriptase inhibitor abacavir. *Lancet* 359: 727-732
  94. Ottinger H. D., Ferencik S., Beelen D. W., Lindemann M., Peceny R., et al. (2003) Hematopoietic stem cell transplantation: contrasting the outcome of transplantations from HLA-identical siblings, partially HLA-mismatched related donors, and HLA-matched unrelated donors. *Blood* 102: 1131-1137
  95. Wang C., Krishnakumar S., Wilhelmy J., Babrzadeh F., Stepanyan L., et al. (2012) High-throughput, high-fidelity HLA genotyping with deep sequencing. *Proc Natl Acad Sci U S A* 109: 8676-8681
-

## References

Wang C., Krishnakumar S., Wilhelmy J., Babrzadeh F., Stepanyan L., et al. (2012) High-throughput, high-fidelity HLA genotyping with deep sequencing. *Proc Natl Acad Sci U S A* 109: 8676-8681

The authors propose a high-throughput HLA genotyping method that employs the use of a single long range PCR to amplify the genomic DNA spanning the majority of the coding regions of four polymorphic HLA genes (HLA-A, -B, -C, and -DRB1). This extensive coverage method enhances allelic resolution, which reduces combination ambiguity that results from off-phase heterozygous sequences. Their coverage of non-polymorphic regions increases the chance of identifying previously undescribed alleles with mismatches, insertions, and deletions.

Illumina Technology: MiSeq and HiSeq 2000 for 150 and 100 bp paired-end sequencing. *Used Genome Analyzer<sub>IIx</sub> to sequence 150 bases from both ends*

Zheng X., Shen J., Cox C., Wakefield J. C., Ehm M. G., et al. (2014) HIBAG-HLA genotype imputation with attribute bagging. *Pharmacogenomics J* 14: 192-200

Hosomichi K., Jinam T. A., Mitsunaga S., Nakaoka H. and Inoue I. (2013) Phase-defined complete sequencing of the HLA genes by next-generation sequencing. *BMC Genomics* 14: 355

Kim Y. J., Kim H. Y., Lee J. H., Yu S. J., Yoon J. H., et al. (2013) A genome-wide association study identified new variants associated with the risk of chronic hepatitis B. *Hum Mol Genet* 22: 4233-4238


Lam T. H., Shen M., Chia J. M., Chan S. H. and Ren E. C. (2013) Population-specific recombination sites within the human MHC region. *Heredity (Edinb)* 111: 131-138

Levin A. M., Iannuzzi M. C., Montgomery C. G., Trudeau S., Datta I., et al. (2013) Association of ANXA11 genetic variation with sarcoidosis in African Americans and European Americans. *Genes Immun* 14: 13-18

Riolobos L., Hirata R. K., Turtle C. J., Wang P. R., Gornalusse G. G., et al. (2013) HLA Engineering of Human Pluripotent Stem Cells. *Mol Ther* 21: 1232-1241

## Phase-Defined HLA Sequencing

Phase-defined sequencing indicates which of the two parental chromosomes a particular allele is derived from. Paired-end sequencing inherently provides information for phasing. When a read encompasses two or more heterozygous genotypes of an individual, the phase of the heterozygous genotypes is determined since each fragment from which a read or pair of reads is obtained in a single allele. Therefore, if read lengths have sufficiently high coverage a substantial amount of phase information can be obtained.<sup>96</sup> This has vast implications in understanding the interplay of genetic variation and disease,<sup>97</sup> imputing untyped genetic variation,<sup>98-100</sup> calling genotypes in sequence data,<sup>101-104</sup> detecting genotype error,<sup>105</sup> inferring human demographic history,<sup>106</sup> inferring points of recombination,<sup>107</sup> detecting recurrent mutation,<sup>107</sup> signatures of selection,<sup>108</sup> and modeling cis-regulation of gene expression.

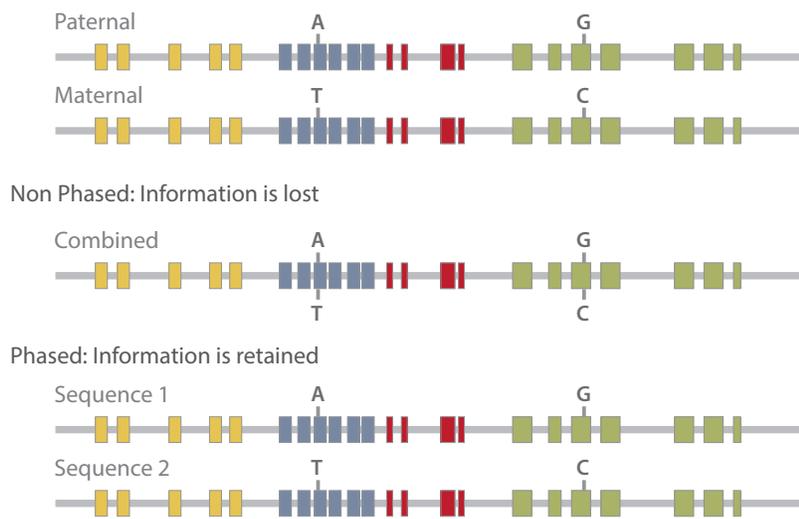

Combination ambiguity occurs when a non-phased consensus sequence is generated. When a single consensus sequence is generated this conceals which of the parental chromosome variants are derived from. Phasing analysis enables researchers to generate two identifiable sequences that correspond to both parental chromosomes. This resolves combination ambiguity and enables researchers to identify which of the two parental chromosomes variants are derived from.

96. Levy S., Sutton G., Ng P. C., Feuk L., Halpern A. L., et al. (2007) The diploid genome sequence of an individual human. *PLoS Biol* 5: e254
97. Tewhey R., Bansal V., Torkamani A., Topol E. J. and Schork N. J. (2011) The importance of phase information for human genomics. *Nat Rev Genet* 12: 215-223
98. Marchini J., Howie B., Myers S., McVean G. and Donnelly P. (2007) A new multipoint method for genome-wide association studies by imputation of genotypes. *Nat Genet* 39: 906-913
99. Browning B. L. and Browning S. R. (2009) A unified approach to genotype imputation and haplotype-phase inference for large data sets of trios and unrelated individuals. *Am J Hum Genet* 84: 210-223
100. Li Y., Willer C. J., Ding J., Scheet P. and Abecasis G. R. (2010) MaCH: using sequence and genotype data to estimate haplotypes and unobserved genotypes. *Genet Epidemiol* 34: 816-834
101. Yu Z., Garner C., Ziogas A., Anton-Culver H. and Schaid D. J. (2009) Genotype determination for polymorphisms in linkage disequilibrium. *BMC Bioinformatics* 10: 63
102. Genomes Project C. (2010) A map of human genome variation from population-scale sequencing. *Nature* 467: 1061-1073
103. Le S. Q. and Durbin R. (2011) SNP detection and genotyping from low-coverage sequencing data on multiple diploid samples. *Genome Res* 21: 952-960
104. Li Y., Sidore C., Kang H. M., Boehnke M. and Abecasis G. R. (2011) Low-coverage sequencing: Implications for design of complex trait association studies. *Genome Res* 21: 940-951
105. Scheet P. and Stephens M. (2008) Linkage disequilibrium-based quality control for large-scale genetic studies. *PLoS Genet* 4: e1000147
106. Tishkoff S. A., Dietsch E., Speed W., Pakstis A. J., Kidd J. R., et al. (1996) Global patterns of linkage disequilibrium at the CD4 locus and modern human origins. *Science* 271: 1380-1387
107. Kong A., Masson G., Frigge M., Gylfason A., Zusmanovich P., et al. (2008) Detection of sharing by descent, long-range phasing and haplotype imputation. *Nat Genet* 40: 1068-1075
108. Sabeti P. C., Reich D. E., Higgins J. M., Levine H. Z., Richter D. J., et al. (2002) Detecting recent positive selection in the human genome from haplotype structure. *Nature* 419: 832-837

Although there has been a strong emergence of next generation sequencing efforts in HLA-genotyping, these comprehensive analyses omit non-coding HLA regions and mRNA-spliced data,<sup>109-111</sup> which may have an impact on gene regulation.<sup>112,113</sup> Moreover, allele determination is conventionally based on sequence alignment to the reference library of HLA sequences in the IMGT/HLA database,<sup>114</sup> which prevents the identification of novel phase-defined HLA gene haplotypes.

## Reviews

Browning S. R. and Browning B. L. (2011) Haplotype phasing: existing methods and new developments. *Nat Rev Genet* 12: 703-714

Tewhey R., Bansal V., Torkamani A., Topol E. J. and Schork N. J. (2011) The importance of phase information for human genomics. *Nat Rev Genet* 12: 215-223

## References

**Hosomichi K., Jinam T. A., Mitsunaga S., Nakaoka H. and Inoue I. (2013) Phase-defined complete sequencing of the HLA genes by next-generation sequencing. *BMC Genomics* 14: 355**

This is the first study to report a complete sequence of the HLA region. Here the authors were able to determine the phase-defined entire HLA gene sequences, regardless of whether the alleles were rare or novel. They sequenced long-range PCR products of HLA genes spanning from promoter to 3'-UTRs and employed a gene-tagging method to generate two HLA gene haplotype sequences based on phase-defined SNVs. Paired end reads of 2 x 250 bps allowed them to demonstrate phase-defined allele determination for 33 HLA homozygous samples, 11 HLA heterozygous samples, and 3 parent-child families.

Illumina Technology: MiSeq 2 x 250 bp and Nextera DNA Sample Prep Kit for library construction


- 
109. Elsner H. A., Bernard G., Eiz-Vesper B., de Matteis M., Bernard A., et al. (2002) Non-expression of HLA-A\*2901102 N is caused by a nucleotide exchange in the mRNA splicing site at the beginning of intron 4. *Tissue Antigens* 59: 139-141
  110. Tamouza R., El Kassab N., Schaeffer V., Carbonnelle E., Tatari Z., et al. (2000) A novel HLA-B\*39 allele (HLA-B\*3916) due to a rare mutation causing cryptic splice site activation. *Hum Immunol* 61: 467-473
  111. Dubois V., Tiercy J. M., Labonne M. P., Dormoy A. and Gebuhrer L. (2004) A new HLA-B44 allele (B\*44020102S) with a splicing mutation leading to a complete deletion of exon 5. *Tissue Antigens* 63: 173-180
  112. Cocco E., Meloni A., Murru M. R., Corongiu D., Tranquilli S., et al. (2012) Vitamin D responsive elements within the HLA-DRB1 promoter region in Sardinian multiple sclerosis associated alleles. *PLoS One* 7: e41678
  113. Thomas R., Apps R., Qi Y., Gao X., Male V., et al. (2009) HLA-C cell surface expression and control of HIV/AIDS correlate with a variant upstream of HLA-C. *Nat Genet* 41: 1290-1294
  114. Hosomichi K., Jinam T. A., Mitsunaga S., Nakaoka H. and Inoue I. (2013) Phase-defined complete sequencing of the HLA genes by next-generation sequencing. *BMC Genomics* 14: 355
-

# SELF VS NON-SELF ANTIGEN DISCRIMINATION

## Tolerance

Tolerance refers to the many layers of protection imposed by the immune system to prevent the reaction of its cells and antibodies against host components. An important form of tolerance is self-tolerance, which refers to the lack of response of the immune system to self-antigens.

Recent studies report a more active role of immune cells in the selective inhibition of responses to self-antigens. For example the study of regulatory T cells (TREG), which in fact recognize self-proteins, have revolutionized the field of tolerance and autoimmunity, not to mention transplantation.

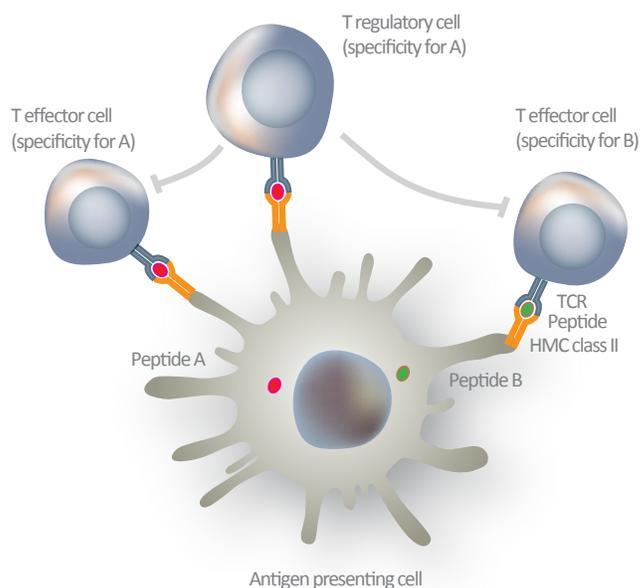

Linked suppression represents a way in which regulatory T cells (TREG) support local self-tolerance. TREG cells inhibit antigen-presenting cells (APCs) presenting their cognate antigen. They can also inhibit bystander T cells, of the same and different antigen specificity, through soluble inhibitory factors.

Central tolerance occurs in the primary lymphoid organs: the bone marrow for B cells and the thymus for T cells. In the first step of this process, T or B-cell clones that recognize self-antigens with high affinity are not allowed to mature. Peripheral tolerance is a secondary precaution in the event that some self-reactive lymphocytes do find their way into the periphery and secondary lymphoid tissues. The peripheral tolerance will render some self-reactive lymphocytes in secondary lymphoid tissues inactive and generates others that actively inhibit immune responses against self. Furthermore, induced cell death, or apoptosis adds a further protective measure by limiting the lifespan of activated lymphocytes.

T cell anergy has been characterized as a hyporesponsive state, or unresponsiveness to an antigenic stimulus, induced by TCR engagement in the absence of costimulation.<sup>115,116</sup> Conversely, when the same antigen is presented with appropriate costimulatory molecules it can become a potent immunogen and mount an immunologic response. Indirect evidence suggests that T cell dysfunction in the tumor microenvironment and establishment of transplant tolerance is partially attributed to T cell anergy.<sup>117</sup> Despite the advances in the characterization of T cell anergy, there are gaps in our knowledge base of the anergic phenotype. This is due to the lack of surface markers that might be useful in identifying anergic T cells. Furthermore, it is unclear teleologically why T cells that are subjected to anergy-inducing conditions are not deleted from the repertoire, in order to eliminate T cells of undesired specificities.

## Reviews

Nishikawa H. and Sakaguchi S. (2014) Regulatory T cells in cancer immunotherapy. *Curr Opin Immunol* 27C: 1-7

Wherry E. J. (2011) T cell exhaustion. *Nat Immunol* 12: 492-499 *Immunol* 132: 170-181

## References

Zheng Y., Zha Y., Spaapen R. M., Mathew R., Barr K., et al. (2013) Egr2-dependent gene expression profiling and ChIP-Seq reveal novel biologic targets in T cell anergy. *Mol Immunol* 55: 283-291

T-cell anergy contributes to peripheral tolerance and plays a role in tumor growth and in promoting transplant allograft acceptance. Egr2 is a critical transcriptional regulator of T-cell anergy. To identify the direct transcriptional targets of Egr2, the authors performed a ChIP-Seq analysis of anti-Egr2 Ab immunoprecipitated nuclear extracts derived from untreated and anergized Th1 T cells. By merging this data with gene expression profiling analyses, the authors revealed 49 targets that are directly regulated by Egr2. They unexpectedly identified cell surface molecules and secreted factors, including lymphocyte-activation gene 3 (Lag3), Class-I-MHC-restricted T-cell associated molecule (Crtam), Semaphorin 7A (Sema7A), and chemokine CCL1. These results indicate that the anergic state may have a functional role through interactions with other immune cells.

Illumina Technology: Genome Analyzer<sub>IIx</sub> for ChIP-Seq of DNA fragments ranging from 200 and 400 bp

Roychoudhuri R., Hirahara K., Mousavi K., Clever D., Klebanoff C. A., et al. (2013) BACH2 represses effector programs to stabilize T(reg)-mediated immune homeostasis. *Nature* 498: 506-510

Genetic polymorphisms within a single locus encoding the transcription factor BACH2 are associated with numerous autoimmune and allergic diseases. The authors use mice that had the BACH2 gene disrupted. They find that BACH2 is a regulator of immune activation.

Illumina Technology: HiSeq for MiRNA sequencing and ChIP-Seq

Ferraro A., D'Alise A. M., Raj T., Asinovski N., Phillips R., et al. (2014) Interindividual variation in human T regulatory cells. *Proc Natl Acad Sci U S A* 111: E1111-1120

- 
115. Steinert E. M., Schwartz R. H. and Singh N. J. (2012) At low precursor frequencies, the T-cell response to chronic self-antigen results in anergy without deletion. *Eur J Immunol* 42: 2875-2880
  116. Chappert P. and Schwartz R. H. (2010) Induction of T cell anergy: integration of environmental cues and infectious tolerance. *Curr Opin Immunol* 22: 552-559
  117. Gajewski T. F., Fuertes M., Spaapen R., Zheng Y. and Kline J. (2011) Molecular profiling to identify relevant immune resistance mechanisms in the tumor microenvironment. *Curr Opin Immunol* 23: 286-292
-

# Autoimmunity

Autoimmunity is caused from the failure of tolerance to protect the host from autoreactive T or B-cell clones. The pathogenesis of these diseases manifests itself in the destruction of proteins, cells, and organs by self-reactive lymphocytes. The onset and pathogenesis of autoimmunity not only depends on intrinsic factors of T and B-cells, such as germline or somatic mutations<sup>118-120</sup> but also on environmental factors such as microbiota or infections,<sup>121</sup> the cytokine milieu, and the presence of other immune cells in the microenvironment.<sup>122</sup>

For most chronic autoimmune and inflammatory diseases, patient populations are heterogeneous and do not uniformly respond to a given therapy. As a result the therapeutic decisions for most autoimmune and inflammatory diseases are based mainly on trial-and-error observations. The development of “actionable biomarkers” may potentially improve the design of clinical trials and inform treatment decisions.<sup>123</sup> For example high-throughput DNA sequencing facilitates the tacking of disease-associated clones of T and B-cells in autoimmune diseases. Furthermore, changes in these cell populations can be correlated with a patient’s response to therapies.<sup>124</sup>

The greatest hope in treating these diseases lie in a greater understanding of the functional roles of genetic and epigenetic variants in autoimmune pathogenesis.<sup>125</sup> Next-generation sequencing of whole-exomes and whole-genomes has become an essential tool in identifying rare genetic variants in large cohorts of autoimmune disease patients. In addition, recent advances in the sequencing of epigenetic markers are adding more information to elucidate the subtle interplay among epigenetic modifications, genetic factors, and environmental signals that predispose individuals to autoimmune risk.

## Reviews

Chan A. C. and Behrens T. W. (2013) Personalizing medicine for autoimmune and inflammatory diseases. *Nat Immunol* 14: 106-109

Graham D. B. and Xavier R. J. (2013) From genetics of inflammatory bowel disease towards mechanistic insights. *Trends Immunol* 34: 371-378

Knoechel B. and Lohr J. G. (2013) Genomics of lymphoid malignancies reveal major activation pathways in lymphocytes. *J Autoimmun* 45: 15-23

Orrù V., Steri M., Sole G., Sidore C., Virdis F., et al. (2013) Genetic variants regulating immune cell levels in health and disease. *Cell* 155: 242-256

Viatte S., Plant D. and Raychaudhuri S. (2013) Genetics and epigenetics of rheumatoid arthritis. *Nat Rev Rheumatol* 9: 141-153

Maecker H. T., Lindstrom T. M., Robinson W. H., Utz P. J., Hale M., et al. (2012) New tools for classification and monitoring of autoimmune diseases. *Nat Rev Rheumatol* 8: 317-328

- 
118. Corradin O., Saiakhova A., Akhtar-Zaidi B., Myeroff L., Willis J., et al. (2014) Combinatorial effects of multiple enhancer variants in linkage disequilibrium dictate levels of gene expression to confer susceptibility to common traits. *Genome Res* 24: 13-Jan
  119. Li Y., Cheng H., Zuo X. B., Sheng Y. J., Zhou F. S., et al. (2013) Association analyses identifying two common susceptibility loci shared by psoriasis and systemic lupus erythematosus in the Chinese Han population. *J Med Genet* 50: 812-818
  120. McDonald-McGinn D. M., Fahiminiya S., Revil T., Nowakowska B. A., Suhl J., et al. (2013) Hemizygous mutations in SNAP29 unmask autosomal recessive conditions and contribute to atypical findings in patients with 22q11.2DS. *J Med Genet* 50: 80-90
  121. Yurkovetskiy L., Burrows M., Khan A. A., Graham L., Volchkov P., et al. (2013) Gender bias in autoimmunity is influenced by microbiota. *Immunity* 39: 400-412
  122. Knoechel B. and Lohr J. G. (2013) Genomics of lymphoid malignancies reveal major activation pathways in lymphocytes. *J Autoimmun* 45: 15-23
  123. Chan A. C. and Behrens T. W. (2013) Personalizing medicine for autoimmune and inflammatory diseases. *Nat Immunol* 14: 106-109
  124. Maecker H. T., Lindstrom T. M., Robinson W. H., Utz P. J., Hale M., et al. (2012) New tools for classification and monitoring of autoimmune diseases. *Nat Rev Rheumatol* 8: 317-328
  125. Laird, P. W. (2010) Principles and challenges of genome-wide DNA methylation analysis. *Nature reviews. Genetics* 11: 191–203
-

## References

**Joseph C. G., Darrach E., Shah A. A., Skora A. D., Casciola-Rosen L. A., et al. (2014) Association of the autoimmune disease scleroderma with an immunologic response to cancer. *Science* 343: 152-157**

Scleroderma is an autoimmune connective tissue disease in which patients make antibodies to a limited group of autoantigens. Patients with scleroderma and antibodies against RPC1 are at increased risk for cancer. The authors sequenced the tumor and normal coding sequences of the POLR3A, TOP1, and CENPB genes in 16 patients. The results suggest that POLR3A mutations triggered cellular immunity and cross-reactive humoral immune responses.

Illumina Technology: Genome Analyzer<sub>IIx</sub>

**Christodoulou K., Wisikin A. E., Gibson J., Tapper W., Willis C., et al. (2013) Next generation exome sequencing of paediatric inflammatory bowel disease patients identifies rare and novel variants in candidate genes. *Gut* 62: 977-984**

The authors utilize exome-sequencing analysis to identify rare and novel variants in known inflammatory bowel disease (IBD) susceptibility genes. Of a panel of 169 known IBD susceptibility genes, approximately 300 non-synonymous, truncating and frameshift mutations were identified from eight pediatric IBD patients. After excluding HLA variants, they uncovered 58 variants across 39 genes, of which 17 were not previously reported. Of the cohort, both patients with severe ulcerative colitis (UC) displayed a distinct profile; both carried potentially deleterious unique variation in the B-cell regulatory gene BACH2 and IL10 genes, which was not seen in the other IBD patients. Variation in BACH2 has not been reported in GWAS of UC.

Illumina Technology: HiSeq 2000 for exome sequencing

**Coit P., Jeffries M., Aitorok N., Dozmorov M. G., Koelsch K. A., et al. (2013) Genome-wide DNA methylation study suggests epigenetic accessibility and transcriptional poisoning of interferon-regulated genes in naive CD4+ T cells from lupus patients. *J Autoimmun* 43: 78-84**

DNA methylation changes have been implicated in T cell differentiation in systemic lupus erythematosus patients. The authors performed a genome-wide DNA methylation study and gene expression profiling in naive CD4+ T cells from lupus patients and controls. Among 86 CG sites that are differentially methylated in naive CD4+ T cells from lupus patients, they revealed that 21 out of 35 hypo methylated genes are regulated by type-1 interferon - including IFIT1, IFIT3, MX1, STAT1, IFI44L, USP18, TRIM22 and BST2. These results indicate that abnormal DNA methylation exists in lupus T cells prior to activation and differentiation and provide an epigenetic explanation for hyper-responsiveness to type-1 interferon in lupus T cells.

Illumina Technology: Infinium Human-Methylation450 BeadChip array for DNA methylation studies and HumanHT-12 v4 Expression BeadChip array for gene expression studies. DNA methylation analysis was performed with GenomeStudio methylation analysis package

**Han A., Newell E. W., Glanville J., Fernandez-Becker N., Khosla C., et al. (2013) Dietary gluten triggers concomitant activation of CD4+ and CD8+ alpha beta T cells and gamma delta T cells in celiac disease. *Proc Natl Acad Sci U S A* 110: 13073-13078**

Celiac disease is an intestinal autoimmune disease caused by dietary gluten and gluten-specific CD4+ T cell responses. Gluten exposure also induces the appearance of activated, gut-homing CD8+  $\alpha\beta$  and  $\gamma\delta$  T cells in peripheral blood. Single-cell T cell receptor sequence analysis indicates that both of these cell populations have highly focused T cell receptor repertoires. Such a focused repertoire usually indicates that the induction is driven by an antigen.

Illumina Technology: MiSeq paired-end sequencing

**Roychoudhuri R., Hirahara K., Mousavi K., Clever D., Klebanoff C. A., et al. (2013) BACH2 represses effector programs to stabilize T(reg)-mediated immune homeostasis. *Nature* 498: 506-510**

BACH2 is expressed in B cells where it acts as a transcriptional repressor of Blimp-1 and other class switch recombination genes. Polymorphisms within a single locus encoding the transcription factor BACH2 are associated with numerous autoimmune and allergic diseases. By studying mice in which the BACH2 gene was disrupted, the authors found that BACH2 is a key regulator of CD4+ T-cell differentiation that prevents inflammatory disease by controlling the balance between tolerance and immunity. They stimulated BACH2 knockout murine naive CD4+ T cells, and subsequently performed massively parallel RNA sequencing to show that the majority of differentially expressed genes were unregulated in BACH2-deficient cells. The authors measured genome-wide BACH2 binding in iTreg cells by chromatin immunoprecipitation with massively parallel sequencing. They determined that BACH2 bound 43.6% of all derepressed genes, including 408 derepressed effector lineage-associated genes.

Illumina Technology: HiSeq 2000 and used TruSeq Sample Prep Kit to prepare RNA-Seq libraries

**Lessard C. J., Li H., Adrianto I., Ice J. A., Rasmussen A., et al. (2013) Variants at multiple loci implicated in both innate and adaptive immune responses are associated with Sjögren's syndrome. Nat Genet 45: 1284-1292**

Sjögren's syndrome is a common autoimmune disease that typically presents as inflammation of the cornea (keratoconjunctivitis sicca) and dry mouth syndrome (xerostomia). In this publication the authors performed a genome-wide association study using Illumina Omni1 Quad array, Illumina ImmunoChip and performed gene expression profiling on the Illumina Human WG-6 v.3.0 BeadChip. Using bioinformatics... read more tools, the combination of genome-wide significance thresholding and suggestively associated variants provided evidence of direct and indirect protein-protein interaction and enrichment of genes involved in immune signaling processes including: TNFAIP3, PTTG1, PRDM1, DGKQ, FCGR2A, IRAK1BP1, ITSN2 and PHIP, among others.

Illumina Technology: HumanOmni1-Quad, Human ImmunoChip

**Martin J. E., Assassi S., Diaz-Gallo L. M., Broen J. C., Simeon C. P., et al. (2013) A systemic sclerosis and systemic lupus erythematosus pan-meta-GWAS reveals new shared susceptibility loci. Hum Mol Genet 22: 4021-4029**

In this meta-analysis of two genome-wide association studies (GWAS) the authors searched for common genetic susceptibility loci for systemic sclerosis (SSc) and lupus erythematosus (SLE). Both diseases are archetypical autoimmune diseases and previous studies have shown that several autoimmune diseases have a common genetic basis. In this study of a total of ~21,000 samples, one new associated locus was identified and two previously described SLE loci were found to be shared with SSc.

Illumina Technology: HumanCNV370-Duo, HumanHap550, Human610-Quad, Human Gene Expression - BeadArray

Bhanusali D. G., Sachdev A., Olson M. A., Gerlach J. A. and Sinha A. A. (2014) PTPN22 profile indicates a novel risk group in Alopecia areata. Hum Immunol 75: 81-87

Leung J. M., Davenport M., Wolff M. J., Wiens K. E., Abidi W. M., et al. (2014) IL-22-producing CD4+ cells are depleted in actively inflamed colitis tissue. Mucosal Immunol 7: 124-133

Furukawa H., Oka S., Matsui T., Hashimoto A., Arinuma Y., et al. (2013) Genome, epigenome and transcriptome analyses of a pair of monozygotic twins discordant for systemic lupus erythematosus. Hum Immunol 74: 170-175

Guo X., Brenner M., Zhang X., Laragione T., Tai S., et al. (2013) Whole-genome sequences of da and f344 rats with different susceptibilities to arthritis, autoimmunity, inflammation and cancer. Genetics 194: 1017-1028

Koelsch K. A., Webb R., Jeffries M., Dozmorov M. G., Frank M. B., et al. (2013) Functional characterization of the MECP2/IRAK1 lupus risk haplotype in human T cells and a human MECP2 transgenic mouse. J Autoimmun 41: 168-174

Nakano K., Whitaker J. W., Boyle D. L., Wang W. and Firestein G. S. (2013) DNA methylome signature in rheumatoid arthritis. Ann Rheum Dis 72: 110-117

Wang C., Ahlford A., Laxman N., Nordmark G., Eloranta M. L., et al. (2013) Contribution of IKBKE and IFIH1 gene variants to SLE susceptibility. Genes Immun 14: 217-222

Ahn J., Gutman D., Saijo S. and Barber G. N. (2012) STING manifests self DNA-dependent inflammatory disease. Proc Natl Acad Sci U S A 109: 19386-19391

Cottrell T. R., Hall J. C., Rosen A. and Casciola-Rosen L. (2012) Identification of novel autoantigens by a triangulation approach. J Immunol Methods 385: 35-44

Labbe C., Boucher G., Foisy S., Alikashani A., Nkwimi H., et al. (2012) Genome-wide expression profiling implicates a MAST3-regulated gene set in colonic mucosal inflammation of ulcerative colitis patients. Inflamm Bowel Dis 18: 1072-1080


## Solid Organ Transplantation

Graft rejection in solid organ transplantation is attributed to histoincompatible tissues. One type of transplanted tissue is an allograft, which is transferred between genetically different members of the same species. Because an allograft is genetically dissimilar to the host and therefore expresses unique antigens, these are often not recognized as self-antigens by the immune system and result in graft rejection.

Tissues that share sufficient antigenic similarity, allowing transfer without immunologic rejection, are said to be histoincompatible, as is the case when the transfer occurs between identical twins. Most transplants are conducted between individuals with a matching ABO blood group and HLA matching. However, even when MHC antigens are identical, the transplanted tissue can be rejected because of differences at various other loci, including the minor histocompatibility locus.

Currently, for heart transplant recipients, the endomyocardial biopsy (EMB) has been employed as the 'gold standard' for rejection surveillance. However, the endomyocardial biopsy is an expensive and invasive procedure that is limited by sampling error, interobserver variability in grading, late detection of rejection, and risk of morbidity.<sup>126,127</sup> Therefore, there has been a considerable effort to develop noninvasive techniques that might replace or reduce the need for EMB, with much focus placed on monitoring the recipient's immune response to detect the onset of rejection.

- 
126. Winters G. L. and McManus B. M. (1996) Consistencies and controversies in the application of the International Society for Heart and Lung Transplantation working formulation for heart transplant biopsy specimens. Rapamycin Cardiac Rejection Treatment Trial Pathologists. *J Heart Lung Transplant* 15: 728-735
127. Oto T., Levvey B. J. and Snell G. I. (2007) Potential refinements of the International Society for Heart and Lung Transplantation primary graft dysfunction grading system. *J Heart Lung Transplant* 26: 431-436
- 

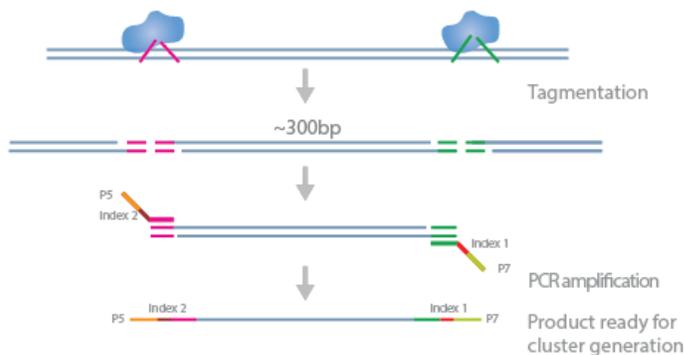

Gene Name  
**HLA-A\*02:10:01:02N**

Subtype  
 Non-synonymous substitutions  
 Non-Coding substitutions

Type  
 Synonymous substitutions in coding region  
 Protein expression

## Reviews

Gabriel C., Furst D., Fae I., Wenda S., Zollikofer C., et al. (2014) HLA typing by next-generation sequencing - getting closer to reality. *Tissue Antigens* 83: 65-75

Boyd S. D. (2013) Diagnostic applications of high-throughput DNA sequencing. *Annu Rev Pathol* 8: 381-410

De Santis D., Dinauer D., Duke J., Erlich H. A., Holcomb C. L., et al. (2013) 16(th) IHIW : review of HLA typing by NGS. *Int J Immunogenet* 40: 72-76

Starzl R., Brandacher G., Lee W. P., Carbonell J., Zhang W., et al. (2013) Review of the early diagnoses and assessment of rejection in vascularized composite allotransplantation. *Clin Dev Immunol* 2013: 402980

## References

[Chen Y., Zhang H., Xiao X., Jia Y., Wu W., et al. \(2013\) Peripheral blood transcriptome sequencing reveals rejection-relevant genes in long-term heart transplantation. \*Int J Cardiol\* 168: 2726-2733](#)

The authors employed transcriptome sequencing of peripheral blood mononuclear cells (PBMCs) derived from 6 quiescent and 6 severe rejection heart transplant recipients. Through digital gene expression (DGE) profiling, a measurement of expression based on the number of reads of the same or similar sequences, they identified a 10-gene PBMC signature capable of distinguishing patients with acute cardiac allograft rejection. Based on a protein-protein interaction network analysis, the authors indicate that CXCR4 and HLA-A are the most informative genes based on a higher degree of control over information flowing to the other 10 genes in the cooperative network.

Illumina Technology: Genome Analyzer<sub>IIx</sub> for 85 bp reads RNA-Seq

[Hosomichi K., Jinam T. A., Mitsunaga S., Nakaoka H. and Inoue I. \(2013\) Phase-defined complete sequencing of the HLA genes by next-generation sequencing. \*BMC Genomics\* 14: 355](#)

The human leukocyte antigen (HLA) region, the 3.8-Mb segment of the human genome at 6p21, has been associated with more than 100 different diseases, mostly autoimmune diseases. Due to the complex nature of HLA genes, there are difficulties in elucidating complete HLA gene sequences especially HLA gene haplotype structures by the conventional sequencing method. This study presents a new method for... read more cost-effective phase-defined complete sequencing of HLA genes using indexed multiplexed samples on Illumina MiSeq. The method was demonstrated on 53 samples showing high resolution for HLA typing.

Illumina Technology: MiSeq, Nextera DNA Sample Prep

[Snyder T. M., Khush K. K., Valentine H. A. and Quake S. R. \(2011\) Universal noninvasive detection of solid organ transplant rejection. \*Proceedings of the National Academy of Sciences of the United States of America\* 108: 6229-6234](#)

Due to increased cell death in the organ during graft rejection, increased donor molecules are expected to be present in the blood at these times. Here the authors genotyped the donor and recipient to establish a unique donor "genetic fingerprint," which was subsequently detected by high-throughput sequencing of the cell-free DNA in peripheral blood of heart transplant recipients. Reads with donor and recipient SNP calls were identified to determine a % Donor DNA. This study establishes a mean value below 1% as indicative of a healthy normal level of donor-derived cell-free DNA. In contrast, during organ rejection the level of donor DNA signal rises to a mean value ranging from 3-4% of the total cell-free DNA.

Illumina Technology: Genome Analyzer<sub>IIx</sub> and Omni1-Quad Beadchip

# INFECTIOUS DISEASES AND VACCINES

New technological advances in T cell isolation and T receptor sequencing have enabled greater understanding of the basic structure of immune T cell repertoires, the diversity of responses within and between individuals, and temporal changes in repertoires and in response to infectious conditions.

## Reviews

Vladimer G. I., Marty-Roix R., Ghosh S., Weng D. and Lien E. (2013) Inflammasomes and host defenses against bacterial infections. *Curr Opin Microbiol* 16: 23-31

## Viral Infections

Viral infections such as HIV, a retrovirus, are able to perturb and alter gene expression through several mechanisms. Studies have profiled the expression of cellular miRNA and some sncRNA post HIV infection using next generation sequencing.<sup>128-130</sup> Emerging studies have focused on the novel mechanisms of gene expression regulation, central to recently discovered players between HIV and the immune system. For example, the human leukocyte antigen (HLA) family of proteins plays a key role in retroviral progression because it is a crucial modulator of the immune response.

Ultimately, understanding how immune cells, such as naïve virus-specific CD8 T cells, influence the type of immune response generated after virus infections is critical to the development of enhanced therapeutic and vaccination strategies to exploit CD8+ T cell-mediated immunity.

## Reviews

Bauersachs S. and Wolf E. (2013) Immune aspects of embryo-maternal cross-talk in the bovine uterus. *J Reprod Immunol* 97: 20-26

Celsi F., Catamo E., Kleiner G., Tricarico P. M., Vuch J., et al. (2013) HLA-G/C, miRNAs, and their role in HIV infection and replication. *Biomed Res Int* 2013: 693643

La Gruta N. L. and Thomas P. G. (2013) Interrogating the relationship between naïve and immune antiviral T cell repertoires. *Curr Opin Virol* 3: 447-451

Lipkin W. I. and Firth C. (2013) Viral surveillance and discovery. *Curr Opin Virol* 3: 199-204

- 
128. Duskova K., Nagilla P., Le H. S., Iyer P., Thalamuthu A., et al. (2013) MicroRNA regulation and its effects on cellular transcriptome in human immunodeficiency virus-1 (HIV-1) infected individuals with distinct viral load and CD4 cell counts. *BMC Infect Dis* 13: 250
  129. Whisnant A. W., Bogerd H. P., Flores O., Ho P., Powers J. G., et al. (2013) In-depth analysis of the interaction of HIV-1 with cellular microRNA biogenesis and effector mechanisms. *MBio* 4: e000193
  130. Chang S. T., Thomas M. J., Sova P., Green R. R., Palermo R. E., et al. (2013) Next-generation sequencing of small RNAs from HIV-infected cells identifies phased microrna expression patterns and candidate novel microRNAs differentially expressed upon infection. *MBio* 4: e00549-00512
-

## References

**Dillon S. M., Lee E. J., Kotter C. V., Austin G. L., Dong Z., et al. (2014) An altered intestinal mucosal microbiome in HIV-1 infection is associated with mucosal and systemic immune activation and endotoxemia. *Mucosal Immunol***

In this paper the authors investigated the impact of HIV-1 infection on the intestinal microbiome and its association with mucosal T-cell and dendritic cell (DC) frequency and activation, as well as with levels of systemic T-cell activation, inflammation, and microbial translocation. They found that HIV-1-related change in the microbiome that was associated with increased mucosal cellular immune activation, microbial translocation, and blood T-cell activation.

Illumina Technology: MiSeq with 250 bp paired-end kit

**O'Connor K. S., Parnell G., Patrick E., Ahlenstiel G., Suppiah V., et al. (2014) Hepatic metallothionein expression in chronic hepatitis C virus infection is IFNL3 genotype-dependent. *Genes Immun* 15: 88-94**

The IFNL3 genotype predicts the clearance of hepatitis C virus (HCV), spontaneously and with interferon (IFN)-based therapy. The authors identified an association between a cluster of ISGs, the metallothioneins (MTs) and IFNL3 genotype. They found that metallothioneins (MTs) were significantly upregulated (in contrast to most other ISGs) in HCV-infected liver biopsies of IFNL3 genotype rs8099917 responders.

Illumina Technology: HiSeq 2000 TruSeq RNA sample preparation and Human HT-12\_V3

**Wang X., Wang H. K., Li Y., Hafner M., Banerjee N. S., et al. (2014) microRNAs are biomarkers of oncogenic human papillomavirus infections. *Proc Natl Acad Sci U S A* 111: 4262-4267**

The authors studied miRNA expression in 158 cervical specimens, including 38 normal, 52 cervical intraepithelial neoplasia (CIN), and 68 cervical cancer (CC) tissues. They found an increase of miR-25, miR-92a, and miR-378 expression with lesion progression but no obvious change of miR-22, miR-29a, and miR-100 among the HPV-infected tissues. An expression ratio  $\geq 1.5$  of miR-25/92a group over miR-22/29a group could serve as a cutoff value to distinguish normal cervix from CIN and from CIN to CC.

Illumina Technology: HiSeq 2000

**Chang S. T., Thomas M. J., Sova P., Green R. R., Palermo R. E., et al. (2013) Next-generation sequencing of small RNAs from HIV-infected cells identifies phased microRNA expression patterns and candidate novel microRNAs differentially expressed upon infection. *MBio* 4: e00549-00512**

The authors investigated the effects of HIV infection on small RNA expression in CD4-expressing T lymphoblastoid cells at 5, 12, and 24 h post infection (hpi). The authors focused on the host response at the level of a single infected cell type and profiled this system over time to detect a phased pattern of microRNA expression. Small RNA-Seq identified 14 differentially expressed microRNA at 5 and 12 hpi; many of which displayed initial suppressed expression followed by rebound later by 24 hpi. They also identified a novel microRNA, an 18-mer encoded in the first intron of the EPB41L2 gene, which was highly expressed in uninfected cell and down regulated by 90% at 24 hpi.

Illumina Technology: Genome Analyzer<sub>IIx</sub> for RNA-Seq of 54 bp reads. Small RNA libraries were prepared with a small RNA version 1.5 sample preparation kit

**Whisnant A. W., Bogerd H. P., Flores O., Ho P., Powers J. G., et al. (2013) In-depth analysis of the interaction of HIV-1 with cellular microRNA biogenesis and effector mechanisms. *MBio* 4: e00019**

The question of how HIV-1 interfaces with cellular microRNA (miRNA) biogenesis and effector mechanisms has been highly controversial. In this paper, the authors used Illumina HiSeq 2000 for deep sequencing of small RNAs in two different infected cell lines and two types of primary human cells to unequivocally demonstrate that HIV-1 does not encode any viral miRNAs.

Illumina Technology: HiSeq 2000 with TruSeq RNA kit deep sequencing of small RNAs and PAR-CLIP to find miRNA binding sites in the HIV-1 genome

**Genolet R., Leignadier J., Osteras M., Farinelli L., Stevenson B. J., et al. (2014) Duality of the murine CD8 compartment. *Proc Natl Acad Sci U S A* 111: E1007-1015**

**Paquette S. G., Banner D., Chi I., Le?n A. J., Xu L., et al. (2014) Pandemic H1N1 influenza A directly induces a robust and acute inflammatory gene signature in primary human bronchial epithelial cells downstream of membrane fusion. *Virology* 448: 91-103**

**Schoggins J. W., MacDuff D. A., Imanaka N., Gainey M. D., Shrestha B., et al. (2014) Pan-viral specificity of IFN-induced genes reveals new roles for cGAS in innate immunity. *Nature* 505: 691-695**

Bolen C. R., Robek M. D., Brodsky L., Schulz V., Lim J. K., et al. (2013) The blood transcriptional signature of chronic hepatitis C virus is consistent with an ongoing interferon-mediated antiviral response. *J Interferon Cytokine Res* 33: 15-23

Duskova K., Nagilla P., Le H. S., Iyer P., Thalamuthu A., et al. (2013) MicroRNA regulation and its effects on cellular transcriptome in human immunodeficiency virus-1 (HIV-1) infected individuals with distinct viral load and CD4 cell counts. *BMC Infect Dis* 13: 250

Grainger J. R., Wohlfert E. A., Fuss I. J., Bouladoux N., Askenase M. H., et al. (2013) Inflammatory monocytes regulate pathologic responses to commensals during acute gastrointestinal infection. *Nat Med* 19: 713-721

Haralambieva I. H., Oberg A. L., Ovsyannikova I. G., Kennedy R. B., Grill D. E., et al. (2013) Genome-wide characterization of transcriptional patterns in high and low antibody responders to rubella vaccination. *PLoS One* 8: e62149

Pociask D. A., Scheller E. V., Mandalapu S., McHugh K. J., Enelow R. I., et al. (2013) IL-22 Is Essential for Lung Epithelial Repair following Influenza Infection. *Am J Pathol* 182: 1286-1296

Swaminathan S., Hu X., Zheng X., Kriga Y., Shetty J., et al. (2013) Interleukin-27 treated human macrophages induce the expression of novel microRNAs which may mediate anti-viral properties. *Biochem Biophys Res Commun* 434: 228-234

Archer J., Weber J., Henry K., Winner D., Gibson R., et al. (2012) Use of Four Next-Generation Sequencing Platforms to Determine HIV-1 Coreceptor Tropism. *PLoS ONE* 7: e49602

## Vaccine Development

The new generation of sequencing technology holds tremendous promise in the areas of systems biology<sup>131</sup> and vaccinomics<sup>132-136</sup> for developing a deeper understanding of the host response to both vaccines and viral infections.

Studies, which report differential gene expression patterns between high and low responders to vaccines, provide insight into the divergent immunoregulatory processes between high and low responders. Further investigation of these loci may lead to important findings regarding the genetic control of immune responses, which can inform the engineering of new vaccine candidates.<sup>132,136</sup>

Human monoclonal antibodies have a high potential to serve as potential therapeutic tools. Until recently, single antibodies capable of neutralizing a broad array of evolving viruses, such as influenza or HIV, were considered extremely rare and nearly impossible to isolate. By employing high-throughput technologies, careful screening processes and clever selection of infected donors, researchers are now able to isolate and characterize these broadly neutralizing antibodies.<sup>137,138</sup> There is now a strong effort to preferentially target the epitopes of these antibodies.

### Reviews

Wilson P. C. and Andrews S. F. (2012) Tools to therapeutically harness the human antibody response. *Nat Rev Immunol* 12: 709-719

Kaur K., Sullivan M. and Wilson P. C. (2011) Targeting B cell responses in universal influenza vaccine design. *Trends Immunol* 32: 524-531

- 
131. Oberg A. L., Kennedy R. B., Li P., Ovsyannikova I. G. and Poland G. A. (2011) Systems biology approaches to new vaccine development. *Curr Opin Immunol* 23: 436-443
  132. Poland G. A., Jacobson R. M. and Ovsyannikova I. G. (2009) Trends affecting the future of vaccine development and delivery: the role of demographics, regulatory science, the anti-vaccine movement, and vaccinomics. *Vaccine* 27: 3240-3244
  133. Poland G. A., Ovsyannikova I. G., Kennedy R. B., Haralambieva I. H. and Jacobson R. M. (2011) Vaccinomics and a new paradigm for the development of preventive vaccines against viral infections. *OMICS* 15: 625-636
  134. Poland G. A., Ovsyannikova I. G. and Jacobson R. M. (2008) Personalized vaccines: the emerging field of vaccinomics. *Expert Opin Biol Ther* 8: 1659-1667
  135. Poland G. A., Ovsyannikova I. G., Jacobson R. M. and Smith D. I. (2007) Heterogeneity in vaccine immune response: the role of immunogenetics and the emerging field of vaccinomics. *Clin Pharmacol Ther* 82: 653-664
  136. Poland G. A., Kennedy R. B. and Ovsyannikova I. G. (2011) Vaccinomics and personalized vaccinology: is science leading us toward a new path of directed vaccine development and discovery? *PLoS Pathog* 7: e1002344
  137. Ekiert D. C., Bhabha G., Elsliger M. A., Friesen R. H., Jongeneelen M., et al. (2009) Antibody recognition of a highly conserved influenza virus epitope. *Science* 324: 246-251
  138. Sui J., Hwang W. C., Perez S., Wei G., Aird D., et al. (2009) Structural and functional bases for broad-spectrum neutralization of avian and human influenza A viruses. *Nat Struct Mol Biol* 16: 265-273
-

## References

Furman D., Jovic V., Kidd B., Shen-Orr S., Price J., et al. (2013) Apoptosis and other immune biomarkers predict influenza vaccine responsiveness. *Mol Syst Biol* 9: 659

There is a known association between pre-existing antibodies and poor vaccination response, which has been attributed to pre-existing flu-specific memory CD4+ T cells that inhibit antigen-presentation by dendritic cells, and subsequently suppress B-cell response. The authors used whole-genome DNA microarrays to assess the baseline levels of immune parameters that correlate with the hemagglutinin inhibition titer response in influenza vaccination. They reported 109 gene modules, sets of coexpressed genes to which the same set of transcription factors binds. Nine variables could predict the antibody response with 84% accuracy. This is the first study to report an association between apoptosis of reactive memory cells and robust antibody response to a vaccine.

Illumina Technology: HumanHT-12v3 Expression BeadChip and GenomeStudio software

Kennedy R. B., Oberg A. L., Ovsyannikova I. G., Haralambieva I. H., Grill D., et al. (2013) Transcriptomic profiles of high and low antibody responders to smallpox vaccine. *Genes Immun* 14: 277-285

Vaccinia virus (VACV) is an immunologically cross-protective virus that used in the smallpox vaccine. The authors used mRNA-Seq transcriptome profiling to identify host and viral gene expression patterns in peripheral blood mononuclear cells (PBMCs) from smallpox vaccine recipients after VACV stimulation. Of the over 1200 genes that exhibited differential gene expression, they identified a number of chemokine, cytokines, interferon and macrophage-associated genes with significant down regulation upon vaccinia infection. Conversely, they identified upregulation in genes encoding histone, IFN $\beta$ , IFN $\gamma$  and heat-shock proteins. Gene set analysis revealed that genes with lowest expression values were expressed 'late' in the viral life cycle, whereas genes classified as 'early' were expressed at significantly higher levels. The patient cohort of high and low vaccinia-specific responders enabled the identification of differential gene regulation patterns between robust humoral immunity and weaker humoral immune responses.

Illumina Technology: Genome Analyzer<sub>IIx</sub> and Single Read Cluster Generation Kit (v2) and 50 cycle Sequencing Kit (v3). cDNA libraries were created using mRNA-Seq 8 sample prep kit

Fulton R. W., d'Offay J. M. and Eberle R. (2013) Bovine herpesvirus-1: Comparison and differentiation of vaccine and field strains based on genomic sequence variation. *Vaccine* 31: 1471-1479

Larimore K., McCormick M. W., Robins H. S. and Greenberg P. D. (2012) Shaping of human germline IgH repertoires revealed by deep sequencing. *J Immunol* 189: 3221-3230

# TECHNIQUES

## miRNA and noncoding RNAs

Only a small fraction of the transcriptome is translated, leaving most of the transcriptional output as ncRNAs, which is classified into two broad categories: small and long RNAs. MicroRNA (miRNAs) is a sub-class of the small noncoding RNA (ncRNA) family, which are small endogenously expressed molecules that regulate the expression of proteins encoded by their mRNA targets. miRNAs have been associated with central roles in growth, development, and immune response *in vivo*.<sup>139-141</sup> They primarily target gene expression at the post-transcriptional,<sup>142,143</sup> level by adjoining to the RNA-induced silencing complex (RISC), which targets the 3'-untranslated region (3'-UTR) of complementary mRNAs and results in the transcript's repression or degradation.<sup>144,145</sup>

Recent studies have shown that miRNAs have unique expression profiles in cells of the innate and adaptive immune systems, CNS, and cancers.<sup>146-149</sup> Furthermore, new evidence implicates a central role of miRNAs in altering mRNA expression in HIV-target cells in response to viral replication.<sup>150</sup> Improvements in high-throughput sequencing technologies, with respect to depth and sensitivity, are enabling researchers to profile known and novel miRNAs, and identify their exact sequence and length, which provides insights on RNA editing processes and mutational events.<sup>151</sup> This allows researchers to decode the networks of non-coding RNA control in the development of the adaptive and innate immune systems and their functional response.

“Wherever the requirement for miRNAs has been tested in the immune system, essential roles have been found”

Ansel et al. 2013

139. Barnes M. R., Deharo S., Grocock R. J., Brown J. R. and Sanseau P. (2007) The micro RNA target paradigm: a fundamental and polymorphic control layer of cellular expression. *Expert Opin Biol Ther* 7: 1387-1399
140. Carthew R. W. and Sontheimer E. J. (2009) Origins and Mechanisms of miRNAs and siRNAs. *Cell* 136: 642-655
141. Chua J. H., Armugam A. and Jeyaseelan K. (2009) MicroRNAs: biogenesis, function and applications. *Curr Opin Mol Ther* 11: 189-1996
142. Barik S. (2005) Silence of the transcripts: RNA interference in medicine. *J Mol Med (Berl)* 83: 764-773
143. Amariglio N. and Rechavi G. (2007) A-to-I RNA editing: a new regulatory mechanism of global gene expression. *Blood Cells Mol Dis* 39: 151-155
144. Guo H., Ingolia N. T., Weissman J. S. and Bartel D. P. (2010) Mammalian microRNAs predominantly act to decrease target mRNA levels. *Nature* 466: 835-840
145. Lewis B. P., Burge C. B. and Bartel D. P. (2005) Conserved seed pairing, often flanked by adenosines, indicates that thousands of human genes are microRNA targets. *Cell* 120: 15-201
146. Gomase V. S. and Parundekar A. N. (2009) microRNA: human disease and development. *Int J Bioinform Res Appl* 5: 479-500
147. Fritz J. H., Girardin S. E. and Philpott D. J. (2006) Innate immune defense through RNA interference. *Sci STKE* 2006: pe27
148. Carissimi C., Fulci V. and Macino G. (2009) MicroRNAs: novel regulators of immunity. *Autoimmun Rev* 8: 520-524
149. O'Connell R. M., Rao D. S., Chaudhuri A. A. and Baltimore D. (2010) Physiological and pathological roles for microRNAs in the immune system. *Nat Rev Immunol* 10: 111-122
150. Duskova K., Nagilla P., Le H. S., Iyer P., Thalamuthu A., et al. (2013) MicroRNA regulation and its effects on cellular transcriptome in human immunodeficiency virus-1 (HIV-1) infected individuals with distinct viral load and CD4 cell counts. *BMC Infect Dis* 13: 250
151. Peng Z., Cheng Y., Tan B. C., Kang L., Tian Z., et al. (2012) Comprehensive analysis of RNA-Seq data reveals extensive RNA editing in a human transcriptome. *Nat Biotechnol* 30: 253-260



Chromatin isolation by RNA purification (ChIRP-Seq) is a protocol to detect the locations on the genome where non-coding RNAs (ncRNAs), such as long non-coding RNAs (lncRNAs), and their proteins are bound.<sup>152</sup> In this method, samples are first crosslinked and sonicated. Biotinylated tiling oligos are hybridized to the RNAs of interest, and the complexes are captured with streptavidin magnetic beads. After treatment with RNase H the DNA is extracted and sequenced. With deep sequencing the lncRNA/protein interaction site can be determined at single-base resolution.<sup>153</sup>

## Reviews

Minton K. (2014) Antiviral immunity: Editing HLA-E expression. *Nat Rev Immunol* 14: 210-211

Ansel K. M. (2013) RNA regulation of the immune system. *Immunol Rev* 253: 5-11

Bronevetsky Y. and Ansel K. M. (2013) Regulation of miRNA biogenesis and turnover in the immune system. *Immunol Rev* 253: 304-316

Pagani M., Rossetti G., Panzeri I., de Candia P., Bonnal R. J., et al. (2013) Role of microRNAs and long-non-coding RNAs in CD4(+) T-cell differentiation. *Immunol Rev* 253: 82-96

Leong J. W., Sullivan R. P. and Fehniger T. A. (2012) Natural Killer Cell Regulation by MicroRNAs in Health and Disease. *J Biomed Biotechnol* 2012: 632329

152. Chu C., Qu K., Zhong F. L., Artandi S. E. and Chang H. Y. (2011) Genomic maps of long noncoding RNA occupancy reveal principles of RNA-chromatin interactions. *Mol Cell* 44: 667-678

153. Li Z., Chao T. C., Chang K. Y., Lin N., Patil V. S., et al. (2014) The long noncoding RNA THRIL regulates TNF $\alpha$  expression through its interaction with hnRNPL. *Proc Natl Acad Sci U S A* 111: 1002-1007

## References

Li Z., Chao T. C., Chang K. Y., Lin N., Patil V. S., et al. (2014) The long noncoding RNA THRIL regulates TNF $\alpha$  expression through its interaction with hnRNPL. *Proc Natl Acad Sci U S A* 111: 1002-1007

The non-protein-coding parts of the mammalian genome encode thousands of large intergenic non-coding RNAs (lincRNAs). To identify lincRNAs associated with activation of the innate immune response; this study applied custom microarrays and Illumina RNA sequencing for THP1 macrophages. A panel of 159 lincRNAs was found to be differentially expressed following innate activation. Further analysis of the RNA-Seq data revealed that linc1992 was required for expression of many immune-response genes, including cytokines and regulators of TNF-alpha expression.

Illumina Technology: HiSeq 2000

Kirigin F. F., Lindstedt K., Sellars M., Ciofani M., Low S. L., et al. (2012) Dynamic microRNA gene transcription and processing during T cell development. *J Immunol* 188: 3257-3267

The authors used next generation sequencing to construct a comprehensive miRNA atlas of T cell development, which reveals the dynamic nature of miRNA gene transcription and processing throughout this developmental pathway starting from murine hematopoietic stem cells to mature CD4 and CD8 thymocytes. They found that many of the miRNAs that were highly expressed in bone marrow progenitor populations were down regulated by 2-3 orders of magnitude in thymocyte populations. They also used ChIP-Seq and polyA RNA enrichment (RNA-Seq) to map the structures of miRNA genes expressed in CD4 T cells.

Illumina Technology: Genome Analyzer<sub>IIx</sub> system, RNA-Seq and ChIP-Seq

Wang P., Gu Y., Zhang Q., Han Y., Hou J., et al. (2012) Identification of resting and type I IFN-activated human NK cell miRNomes reveals microRNA-378 and microRNA-30e as negative regulators of NK cell cytotoxicity. *J Immunol* 189: 211-221

Authors used next-generation sequencing to perform smRNA expression profiling of human CD56+CD3-Natural Killer (NK) cells during the process of cytokine activation. Of the >200 novel miRNAs identified, they report two abundant miRNA, miRNA-378 and miRNA-30e, are down-regulated in IFN- $\alpha$  activated NK cells. They also found that these two miRNA directly target granzyme B and perforin, respectively, which demonstrates that miR-378 and miR-30e suppress human NK cell cytotoxicity during NK cell activation.


Illumina Technology: smRNA-Seq between 18 and 30 bp

- Bryniarski K., Ptak W., Jayakumar A., Pullmann K., Caplan M. J., et al. (2013) Antigen-specific, antibody-coated, exosome-like nanovesicles deliver suppressor T-cell microRNA-150 to effector T cells to inhibit contact sensitivity. *J Allergy Clin Immunol* 132: 170-181
- Chang S. T., Thomas M. J., Sova P., Green R. R., Palermo R. E., et al. (2013) Next-generation sequencing of small RNAs from HIV-infected cells identifies phased microRNA expression patterns and candidate novel microRNAs differentially expressed upon infection. *MBio* 4: e00549-00512
- Duskova K., Nagilla P., Le H. S., Iyer P., Thalamuthu A., et al. (2013) MicroRNA regulation and its effects on cellular transcriptome in human immunodeficiency virus-1 (HIV-1) infected individuals with distinct viral load and CD4 cell counts. *BMC Infect Dis* 13: 250
- Pinto S., Michel C., Schmidt-Glenewinkel H., Harder N., Rohr K., et al. (2013) Overlapping gene coexpression patterns in human medullary thymic epithelial cells generate self-antigen diversity. *Proc Natl Acad Sci U S A* 110: E3497-3505
- Su X., Qian C., Zhang Q., Hou J., Gu Y., et al. (2013) miRNomes of haematopoietic stem cells and dendritic cells identify miR-30b as a regulator of Notch1. *Nat Commun* 4: 2903
- Swaminathan S., Hu X., Zheng X., Kriga Y., Shetty J., et al. (2013) Interleukin-27 treated human macrophages induce the expression of novel microRNAs which may mediate anti-viral properties. *Biochem Biophys Res Commun* 434: 228-234
- Whisnant A. W., Bogerd H. P., Flores O., Ho P., Powers J. G., et al. (2013) In-depth analysis of the interaction of HIV-1 with cellular microRNA biogenesis and effector mechanisms. *MBio* 4: e000193
- Liu X., Wang Y., Sun Q., Yan J., Huang J., et al. (2012) Identification of microRNA transcriptome involved in human natural killer cell activation. *Immunol Lett* 143: 208-217
- Nolte-t Hoen E. N., Buermans H. P., Waasdorp M., Stoorvogel W., Wauben M. H., et al. (2012) Deep sequencing of RNA from immune cell-derived vesicles uncovers the selective incorporation of small non-coding RNA biotypes with potential regulatory functions. *Nucleic Acids Res* 40: 9272-9285
- Yang G., Yang L., Zhao Z., Wang J. and Zhang X. (2012) Signature miRNAs Involved in the Innate Immunity of Invertebrates. *PLoS ONE* 7: e39015

## ChIP-Seq

Many transcription factors and chromatin modifiers linked to innate and adaptive immunity.<sup>154-157</sup> The identification and characterization of the genome-wide locations of transcription factors and chromatin-modifying enzymes and the modification status of histones has been accelerated by the application of chromatin immunoprecipitation techniques in next-generation sequencing analysis (ChIP-Seq). This method employs antibodies directed against a target protein to isolate a DNA-protein complex. Purified DNA is obtained from the immunoprecipitated DNA-protein complexes and is subsequently ligated with sequencing adaptors, amplified by PCR and sequenced on a next-generation sequencing platform.<sup>158</sup> Ultimately, the need to comprehend global transcriptional regulation of the immune system positions ChIP-Seq as a powerful application, which informs our understanding of the dynamic processes of stem cell differentiation, formation of immunological memory, disease progression, and response to environmental stimuli.<sup>158,159</sup>

- 
154. Barski A, Cuddapah S, Cui K, Roh T. Y., Schones D. E., et al. (2007) High-resolution profiling of histone methylations in the human genome. *Cell* 129: 823-837
  155. Wei G, Wei L, Zhu J, Zang C, Hu-Li J, et al. (2009) Global mapping of H3K4me3 and H3K27me3 reveals specificity and plasticity in lineage fate determination of differentiating CD4+ T cells. *Immunity* 30: 155-167
  156. Northrup D. L. and Zhao K. (2011) Application of ChIP-Seq and related techniques to the study of immune function. *Immunity* 34: 830-842
  157. Natoli G. (2010) Maintaining cell identity through global control of genomic organization. *Immunity* 33: 12-24
  158. Kidder B. L., Hu G. and Zhao K. (2011) ChIP-Seq: technical considerations for obtaining high-quality data. *Nat Immunol* 12: 918-922
  159. Vahedi G., Takahashi H., Nakayamada S., Sun H. W., Sartorelli V., et al. (2012) STATs Shape the Active Enhancer Landscape of T Cell Populations. *Cell* 151: 981-993
- 

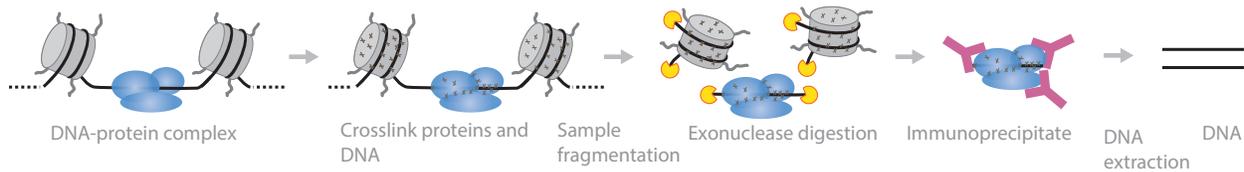

**The ChIP-Seq workflow.** Chromatin immunoprecipitation sequencing (ChIP-Seq) is a well-established method to map specific protein-binding sites. In this method, DNA-protein complexes are crosslinked *in vivo*. Samples are then fragmented and treated with an exonuclease to trim unbound oligonucleotides. Protein-specific antibodies are used to immunoprecipitate the DNA-protein complex. The DNA is extracted and sequenced, giving high-resolution sequences of the protein-binding sites.

## Reviews

Kidder B. L., Hu G. and Zhao K. (2011) ChIP-Seq: technical considerations for obtaining high-quality data. *Nat Immunol* 12: 918-922

## References

Roychoudhuri R., Hirahara K., Mousavi K., Clever D., Klebanoff C. A., et al. (2013) BACH2 represses effector programs to stabilize T(reg)-mediated immune homeostasis. *Nature* 498: 506-510

Wei Y., Chen R., Dimicoli S., Bueso-Ramos C., Neuberger D., et al. (2013) Global H3K4me3 genome mapping reveals alterations of innate immunity signaling and overexpression of JMJD3 in human myelodysplastic syndrome CD34+ cells. *Leukemia* 27: 2177-2186

Zheng Y., Zha Y., Spaapen R. M., Mathew R., Barr K., et al. (2013) Egr2-dependent gene expression profiling and ChIP-Seq reveal novel biologic targets in T cell anergy. *Mol Immunol* 55: 283-291

Hakim O., Resch W., Yamane A., Klein I., Kieffer-Kwon K. R., et al. (2012) DNA damage defines sites of recurrent chromosomal translocations in B lymphocytes. *Nature* 484: 69-74

Kirigin F. F., Lindstedt K., Sellars M., Ciofani M., Low S. L., et al. (2012) Dynamic microRNA gene transcription and processing during T cell development. *J Immunol* 188: 3257-3267

Vahedi G., Takahashi H., Nakayama S., Sun H. W., Sartorelli V., et al. (2012) STATs Shape the Active Enhancer Landscape of T Cell Populations. *Cell* 151: 981-993

Wang P., Gu Y., Zhang Q., Han Y., Hou J., et al. (2012) Identification of resting and type I IFN-activated human NK cell miRNomes reveals microRNA-378 and microRNA-30e as negative regulators of NK cell cytotoxicity. *J Immunol* 189: 211-221

Zhang J. A., Mortazavi A., Williams B. A., Wold B. J. and Rothenberg E. V. (2012) Dynamic Transformations of Genome-wide Epigenetic Marking and Transcriptional Control Establish T Cell Identity. *Cell* 149: 467-482

# BIBLIOGRAPHY

- Ahmed R. and Gray D. (1996) Immunological memory and protective immunity: understanding their relation. *Science* 272: 54-60
- Ahn J., Gutman D., Saijo S. and Barber G. N. (2012) STING manifests self DNA-dependent inflammatory disease. *Proc Natl Acad Sci U S A* 109: 19386-19391
- Amariglio N. and Rechavi G. (2007) A-to-I RNA editing: a new regulatory mechanism of global gene expression. *Blood Cells Mol Dis* 39: 151-155
- Ansel K. M. (2013) RNA regulation of the immune system. *Immunol Rev* 253: 5-11
- Archer J., Weber J., Henry K., Winner D., Gibson R., et al. (2012) Use of Four Next-Generation Sequencing Platforms to Determine HIV-1 Coreceptor Tropism. *PLoS ONE* 7: e49602
- Arthur J. C., Perez-Chanona E., Mühlbauer M., Tomkovich S., Uronis J. M., et al. (2012) Intestinal inflammation targets cancer-inducing activity of the microbiota. *Science* 338: 120-123
- Bae J., Ricciardi C. J., Esposito D., Komarnytsky S., Hu P., et al. (2014) Activation of Pattern Recognition Receptors in Brown Adipocytes Induce Inflammation and Suppress Uncoupling Protein 1 Expression and Mitochondrial Respiration. *Am J Physiol Cell Physiol*
- Barik S. (2005) Silence of the transcripts: RNA interference in medicine. *J Mol Med (Berl)* 83: 764-773
- Barnes M. R., Deharo S., Grocock R. J., Brown J. R. and Sanseau P. (2007) The micro RNA target paradigm: a fundamental and polymorphic control layer of cellular expression. *Expert Opin Biol Ther* 7: 1387-1399
- Barski A., Cuddapah S., Cui K., Roh T. Y., Schones D. E., et al. (2007) High-resolution profiling of histone methylations in the human genome. *Cell* 129: 823-837
- Bauersachs S. and Wolf E. (2013) Immune aspects of embryo-maternal cross-talk in the bovine uterus. *J Reprod Immunol* 97: 20-26
- Benichou J., Ben-Hamo R., Louzoun Y. and Efroni S. (2012) Rep-Seq: uncovering the immunological repertoire through next-generation sequencing. *Immunology* 135: 183-191
- Bhanusali D. G., Sachdev A., Olson M. A., Gerlach J. A. and Sinha A. A. (2014) PTPN22 profile indicates a novel risk group in Alopecia areata. *Hum Immunol* 75: 81-87
- Bolen C. R., Robek M. D., Brodsky L., Schulz V., Lim J. K., et al. (2013) The blood transcriptional signature of chronic hepatitis C virus is consistent with an ongoing interferon-mediated antiviral response. *J Interferon Cytokine Res* 33: 15-23
- Bolotin D. A., Mamedov I. Z., Britanova O. V., Zvyagin I. V., Shagin D., et al. (2012) Next generation sequencing for TCR repertoire profiling: platform-specific features and correction algorithms. *Eur J Immunol* 42: 3073-3083
- Bottazzi B., Doni A., Garlanda C. and Mantovani A. (2010) An integrated view of humoral innate immunity: pentraxins as a paradigm. *Annu Rev Immunol* 28: 157-183
- Bowley D. R., Jones T. M., Burton D. R. and Lerner R. A. (2009) Libraries against libraries for combinatorial selection of replicating antigen-antibody pairs. *Proc Natl Acad Sci U S A* 106: 1380-1385
- Boyd S. D. (2013) Diagnostic applications of high-throughput DNA sequencing. *Annu Rev Pathol* 8: 381-410
- Bronevetsky Y. and Ansel K. M. (2013) Regulation of miRNA biogenesis and turnover in the immune system. *Immunol Rev* 253: 304-316
- Brown E. M., Sadarangani M. and Finlay B. B. (2013) The role of the immune system in governing host-microbe interactions in the intestine. *Nat Immunol* 14: 660-667
- Browning B. L. and Browning S. R. (2009) A unified approach to genotype imputation and haplotype-phase inference for large data sets of trios and unrelated individuals. *Am J Hum Genet* 84: 210-223
- Browning S. R. and Browning B. L. (2011) Haplotype phasing: existing methods and new developments. *Nat Rev Genet* 12: 703-714
- Bryniarski K., Ptak W., Jayakumar A., Pullmann K., Caplan M. J., et al. (2013) Antigen-specific, antibody-coated, exosome-like nanovesicles deliver suppressor T-cell microRNA-150 to effector T cells to inhibit contact sensitivity. *J Allergy Clin Immunol* 132: 170-181
- Caporaso J. G., Lauber C. L., Walters W. A., Berg-Lyons D., Huntley J., et al. (2012) Ultra-high-throughput microbial community analysis on the Illumina HiSeq and MiSeq platforms. *ISME J* 6: 1621-1624
- Carissimi C., Fulci V. and Macino G. (2009) MicroRNAs: novel regulators of immunity. *Autoimmun Rev* 8: 520-524
- Carthew R. W. and Sontheimer E. J. (2009) Origins and Mechanisms of miRNAs and siRNAs. *Cell* 136: 642-655
- Castelino M., Eyre S., Upton M., Ho P. and Barton A. (2013) The bacterial skin microbiome in psoriatic arthritis, an unexplored link in pathogenesis: challenges and opportunities offered by recent technological advances. *Rheumatology (Oxford)*
- Celsi F., Catamo E., Kleiner G., Tricarico P. M., Vuch J., et al. (2013) HLA-G/C, miRNAs, and their role in HIV infection and replication. *Biomed Res Int* 2013: 693643
- Chakalova L. and Fraser P. (2010) Organization of transcription. *Cold Spring Harb Perspect Biol* 2: a000729
- Chan A. C. and Behrens T. W. (2013) Personalizing medicine for autoimmune and inflammatory diseases. *Nat Immunol* 14: 106-109
- Chang S. T., Thomas M. J., Sova P., Green R. R., Palermo R. E., et al. (2013) Next-generation sequencing of small RNAs from HIV-infected cells identifies phased microRNA expression patterns and candidate novel microRNAs differentially expressed upon infection. *MBio* 4: e00549-00512
- Chappert P. and Schwartz R. H. (2010) Induction of T cell anergy: integration of environmental cues and infectious tolerance. *Curr Opin Immunol* 22: 552-559
- Chaumeil J., Micsinai M., Ntziachristos P., Roth D. B., Aifantis I., et al. (2013) The RAG2 C-terminus and ATM protect genome integrity by controlling antigen receptor gene cleavage. *Nat Commun* 4: 2231
- Chen Y., Zhang H., Xiao X., Jia Y., Wu W., et al. (2013) Peripheral blood transcriptome sequencing reveals rejection-relevant genes in long-term heart transplantation. *Int J Cardiol* 168: 2726-2733
- Cho P., Gelinas L., Corbett N. P., Tebbutt S. J., Turvey S. E., et al. (2013) Association of common single-nucleotide polymorphisms in innate immune genes with differences in TLR-induced cytokine production in neonates. *Genes Immun* 14: 199-211
- Christodoulou K., Wiskin A. E., Gibson J., Tapper W., Willis C., et al. (2013) Next generation exome sequencing of paediatric inflammatory bowel disease patients identifies rare and novel variants in candidate genes. *Gut* 62: 977-984
- Chu C., Qu K., Zhong F. L., Artandi S. E. and Chang H. Y. (2011) Genomic maps of long noncoding RNA occupancy reveal principles of RNA-chromatin interactions. *Mol Cell* 44: 667-678

- Chua J. H., Armugam A. and Jeyaseelan K. (2009) MicroRNAs: biogenesis, function and applications. *Curr Opin Mol Ther* 11: 189-199
- Cocco E., Meloni A., Murru M. R., Corongiu D., Tranquilli S., et al. (2012) Vitamin D responsive elements within the HLA-DRB1 promoter region in Sardinian multiple sclerosis associated alleles. *PLoS One* 7: e41678
- Coit P., Jeffries M., Altorok N., Dozmorov M. G., Koelsch K. A., et al. (2013) Genome-wide DNA methylation study suggests epigenetic accessibility and transcriptional poising of interferon-regulated genes in naive CD4+ T cells from lupus patients. *J Autoimmun* 43: 78-84
- Corradin O., Saiakhova A., Akhtar-Zaidi B., Myeroff L., Willis J., et al. (2014) Combinatorial effects of multiple enhancer variants in linkage disequilibrium dictate levels of gene expression to confer susceptibility to common traits. *Genome Res* 24: 13-Jan
- Cottrell T. R., Hall J. C., Rosen A. and Casciola-Rosen L. (2012) Identification of novel autoantigens by a triangulation approach. *J Immunol Methods* 385: 35-44
- Cremer T. and Cremer M. (2010) Chromosome territories. *Cold Spring Harb Perspect Biol* 2: a003889
- Davies J. L., Kawaguchi Y., Bennett S. T., Copeman J. B., Cordell H. J., et al. (1994) A genome-wide search for human type 1 diabetes susceptibility genes. *Nature* 371: 130-136
- de Bakker P. I. and Raychaudhuri S. (2012) Interrogating the major histocompatibility complex with high-throughput genomics. *Hum Mol Genet* 21: R29-36
- De Santis D., Dinauer D., Duke J., Erlich H. A., Holcomb C. L., et al. (2013) 16(th) IHIW : review of HLA typing by NGS. *Int J Immunogenet* 40: 72-76
- Dekosky B. J., Ippolito G. C., Deschner R. P., Lavinder J. J., Wine Y., et al. (2013) High-throughput sequencing of the paired human immunoglobulin heavy and light chain repertoire. *Nat Biotechnol* 31: 166-169
- Dillon S. M., Lee E. J., Kotter C. V., Austin G. L., Dong Z., et al. (2014) An altered intestinal mucosal microbiome in HIV-1 infection is associated with mucosal and systemic immune activation and endotoxemia. *Mucosal Immunol*
- Dubois V., Tiercy J. M., Labonne M. P., Dormoy A. and Gebuhrer L. (2004) A new HLA-B44 allele (B\*44020102S) with a splicing mutation leading to a complete deletion of exon 5. *Tissue Antigens* 63: 173-180
- Duskova K., Nagilla P., Le H. S., Iyer P., Thalamuthu A., et al. (2013) MicroRNA regulation and its effects on cellular transcriptome in human immunodeficiency virus-1 (HIV-1) infected individuals with distinct viral load and CD4 cell counts. *BMC Infect Dis* 13: 250
- Eapen M., Rubinstein P., Zhang M. J., Stevens C., Kurtzberg J., et al. (2007) Outcomes of transplantation of unrelated donor umbilical cord blood and bone marrow in children with acute leukaemia: a comparison study. *Lancet* 369: 1947-1954
- Ekiert D. C., Bhabha G., Elsliger M. A., Friesen R. H., Jongeneelen M., et al. (2009) Antibody recognition of a highly conserved influenza virus epitope. *Science* 324: 246-251
- Elsner H. A., Bernard G., Eiz-Vesper B., de Matteis M., Bernard A., et al. (2002) Non-expression of HLA-A\*2901102 N is caused by a nucleotide exchange in the mRNA splicing site at the beginning of intron 4. *Tissue Antigens* 59: 139-141
- Emerson R., Sherwood A., Desmarais C., Malhotra S., Phippard D., et al. (2013) Estimating the ratio of CD4+ to CD8+ T cells using high-throughput sequence data. *J Immunol Methods* 391: 14-21
- Erlich H. (2012) HLA DNA typing: past, present, and future. *Tissue Antigens* 80: 1-11
- Ferraro A., D'Alise A. M., Raj T., Asinovski N., Phillips R., et al. (2014) Interindividual variation in human T regulatory cells. *Proc Natl Acad Sci U S A* 111: E1111-1120
- Finn J. A. and Crowe J. E., Jr. (2013) Impact of new sequencing technologies on studies of the human B cell repertoire. *Curr Opin Immunol* 25: 613-618
- Fischer N. (2011) Sequencing antibody repertoires: the next generation. *MAbs* 3: 17-20
- Fozza C. and Longinotti M. (2013) T-cell receptor repertoire usage in hematologic malignancies. *Crit Rev Oncol Hematol* 86: 201-211
- Franz B., May K. F., Jr., Dranoff G. and Wucherpfennig K. (2011) Ex vivo characterization and isolation of rare memory B cells with antigen tetramers. *Blood* 118: 348-357
- Fritz J. H., Girardin S. E. and Philpott D. J. (2006) Innate immune defense through RNA interference. *Sci STKE* 2006: pe27
- Fukata M. and Arditi M. (2013) The role of pattern recognition receptors in intestinal inflammation. *Mucosal Immunol* 6: 451-463
- Fulton R. W., d'Offay J. M. and Eberle R. (2013) Bovine herpesvirus-1: Comparison and differentiation of vaccine and field strains based on genomic sequence variation. *Vaccine* 31: 1471-1479
- Furman D., Jovic V., Kidd B., Shen-Orr S., Price J., et al. (2013) Apoptosis and other immune biomarkers predict influenza vaccine responsiveness. *Mol Syst Biol* 9: 659
- Furukawa H., Oka S., Matsui T., Hashimoto A., Arinuma Y., et al. (2013) Genome, epigenome and transcriptome analyses of a pair of monozygotic twins discordant for systemic lupus erythematosus. *Hum Immunol* 74: 170-175
- Gabriel C., Furst D., Fae I., Wenda S., Zollikofer C., et al. (2014) HLA typing by next-generation sequencing - getting closer to reality. *Tissue Antigens* 83: 65-75
- Gajewski T. F., Fustes M., Spaapen R., Zheng Y. and Kline J. (2011) Molecular profiling to identify relevant immune resistance mechanisms in the tumor microenvironment. *Curr Opin Immunol* 23: 286-292
- Genolet R., Stevenson B. J., Farinelli L., Osteras M. and Luescher I. F. (2012) Highly diverse TCRalpha chain repertoire of pre-immune CD8(+) T cells reveals new insights in gene recombination. *EMBO J* 31: 1666-1678
- Genolet R., Leignadier J., Osteras M., Farinelli L., Stevenson B. J., et al. (2014) Duality of the murine CD8 compartment. *Proc Natl Acad Sci U S A* 111: E1007-1015
- Genomes Project C. (2010) A map of human genome variation from population-scale sequencing. *Nature* 467: 1061-1073
- Georgiou G., Ippolito G. C., Beausang J., Busse C. E., Wardemann H., et al. (2014) The promise and challenge of high-throughput sequencing of the antibody repertoire. *Nat Biotechnol* 32: 158-168
- Gomase V. S. and Parundekar A. N. (2009) microRNA: human disease and development. *Int J Bioinform Res Appl* 5: 479-500
- Gottgens B. (2012) Genome-scale technology driven advances to research into normal and malignant haematopoiesis. *Scientifica (Cairo)* 2012: 437956
- Graham D. B. and Xavier R. J. (2013) From genetics of inflammatory bowel disease towards mechanistic insights. *Trends Immunol* 34: 371-378
- Grainger J. R., Wohlfert E. A., Fuss I. J., Bouladoux N., Askenase M. H., et al. (2013) Inflammatory monocytes regulate pathologic responses to commensals during acute gastrointestinal infection. *Nat Med* 19: 713-721
- Guo H., Ingolia N. T., Weissman J. S. and Bartel D. P. (2010) Mammalian microRNAs predominantly act to decrease target mRNA levels. *Nature* 466: 835-840
- Guo X., Brenner M., Zhang X., Laragione T., Tai S., et al. (2013) Whole-genome sequences of da and f344 rats with different susceptibilities to arthritis, autoimmunity, inflammation and cancer. *Genetics* 194: 1017-1028

- Haen S. P. and Rammensee H. G. (2013) The repertoire of human tumor-associated epitopes—identification and selection of antigens and their application in clinical trials. *Curr Opin Immunol* 25: 277-283
- Hakim O., Sung M. H. and Hager G. L. (2010) 3D shortcuts to gene regulation. *Curr Opin Cell Biol* 22: 305-313
- Hakim O., Resch W., Yamane A., Klein I., Kieffer-Kwon K. R., et al. (2012) DNA damage defines sites of recurrent chromosomal translocations in B lymphocytes. *Nature* 484: 69-74
- Han A., Newell E. W., Glanville J., Fernandez-Becker N., Khosla C., et al. (2013) Dietary gluten triggers concomitant activation of CD4+ and CD8+ alphabeta T cells and gammadelta T cells in celiac disease. *Proc Natl Acad Sci U S A* 110: 13073-13078
- Haralambieva I. H., Oberg A. L., Ovsyannikova I. G., Kennedy R. B., Grill D. E., et al. (2013) Genome-wide characterization of transcriptional patterns in high and low antibody responders to rubella vaccination. *PLoS One* 8: e62149
- Holwerda S. J., van de Werken H. J., Ribeiro de Almeida C., Bergen I. M., de Bruijn M. J., et al. (2013) Allelic exclusion of the immunoglobulin heavy chain locus is independent of its nuclear localization in mature B cells. *Nucleic Acids Res* 41: 6905-6916
- Hoogenboom H. R., Griffiths A. D., Johnson K. S., Chiswell D. J., Hudson P., et al. (1991) Multi-subunit proteins on the surface of filamentous phage: methodologies for displaying antibody (Fab) heavy and light chains. *Nucleic Acids Res* 19: 4133-4137
- Hooper L. V., Littman D. R. and Macpherson A. J. (2012) Interactions between the microbiota and the immune system. *Science* 336: 1268-1273
- Horton R., Wilming L., Rand V., Lovering R. C., Bruford E. A., et al. (2004) Gene map of the extended human MHC. *Nat Rev Genet* 5: 889-899
- Hosomichi K., Jinam T. A., Mitsunaga S., Nakaoka H. and Inoue I. (2013) Phase-defined complete sequencing of the HLA genes by next-generation sequencing. *BMC Genomics* 14: 355
- Jankovic M., Feldhahn N., Oliveira T. Y., Silva I. T., Kieffer-Kwon K. R., et al. (2013) 53BP1 alters the landscape of DNA rearrangements and suppresses AID-induced B cell lymphoma. *Mol Cell* 49: 623-631
- Joseph C. G., Darrah E., Shah A. A., Skora A. D., Casciola-Rosen L. A., et al. (2014) Association of the autoimmune disease scleroderma with an immunologic response to cancer. *Science* 343: 152-157
- Kamada N., Chen G. Y., Inohara N. and Nunez G. (2013) Control of pathogens and pathobionts by the gut microbiota. *Nat Immunol* 14: 685-690
- Kamada N., Seo S. U., Chen G. Y. and Nunez G. (2013) Role of the gut microbiota in immunity and inflammatory disease. *Nat Rev Immunol* 13: 321-335
- Kaur K., Sullivan M. and Wilson P. C. (2011) Targeting B cell responses in universal influenza vaccine design. *Trends Immunol* 32: 524-531
- Kawai T. and Akira S. (2011) Toll-like receptors and their crosstalk with other innate receptors in infection and immunity. *Immunity* 34: 637-650
- Kemp C., Mueller S., Goto A., Barbier V., Paro S., et al. (2013) Broad RNA interference-mediated antiviral immunity and virus-specific inducible responses in *Drosophila*. *J Immunol* 190: 650-658
- Kennedy R. B., Oberg A. L., Ovsyannikova I. G., Haralambieva I. H., Grill D., et al. (2013) Transcriptomic profiles of high and low antibody responders to smallpox vaccine. *Genes Immun* 14: 277-285
- Kidder B. L., Hu G. and Zhao K. (2011) ChIP-Seq: technical considerations for obtaining high-quality data. *Nat Immunol* 12: 918-922
- Kim Y. J., Kim H. Y., Lee J. H., Yu S. J., Yoon J. H., et al. (2013) A genome-wide association study identified new variants associated with the risk of chronic hepatitis B. *Hum Mol Genet* 22: 4233-4238
- Kirigin F. F., Lindstedt K., Sellars M., Ciofani M., Low S. L., et al. (2012) Dynamic microRNA gene transcription and processing during T cell development. *J Immunol* 188: 3257-3267
- Klein I. A., Resch W., Jankovic M., Oliveira T., Yamane A., et al. (2011) Translocation-capture sequencing reveals the extent and nature of chromosomal rearrangements in B lymphocytes. *Cell* 147: 95-106
- Knoechel B. and Lohr J. G. (2013) Genomics of lymphoid malignancies reveal major activation pathways in lymphocytes. *J Autoimmun* 45: 15-23
- Koelsch K. A., Webb R., Jeffries M., Dozmorov M. G., Frank M. B., et al. (2013) Functional characterization of the MECP2/IRAK1 lupus risk haplotype in human T cells and a human MECP2 transgenic mouse. *J Autoimmun* 41: 168-174
- Kong A., Masson G., Frigge M., Gylfason A., Zusmanovich P., et al. (2008) Detection of sharing by descent, long-range phasing and haplotype imputation. *Nat Genet* 40: 1068-1075
- Kovalchuk A. L., Ansarah-Sobrinho C., Hakim O., Resch W., Tolarova H., et al. (2012) Mouse model of endemic Burkitt translocations reveals the long-range boundaries of Ig-mediated oncogene deregulation. *Proc Natl Acad Sci U S A* 109: 10972-10977
- La Gruta N. L. and Thomas P. G. (2013) Interrogating the relationship between naive and immune antiviral T cell repertoires. *Curr Opin Virol* 3: 447-451
- Labbe C., Boucher G., Foisy S., Alikashani A., Nkwimi H., et al. (2012) Genome-wide expression profiling implicates a MAST3-regulated gene set in colonic mucosal inflammation of ulcerative colitis patients. *Inflamm Bowel Dis* 18: 1072-1080
- Lagerkvist A. C., Furebring C. and Borrebaeck C. A. (1995) Single, antigen-specific B cells used to generate Fab fragments using CD40-mediated amplification or direct PCR cloning. *Biotechniques* 18: 862-869
- Lam T. H., Shen M., Chia J. M., Chan S. H. and Ren E. C. (2013) Population-specific recombination sites within the human MHC region. *Heredity (Edinb)* 111: 131-138
- Lanier L. L. (2013) Shades of grey--the blurring view of innate and adaptive immunity. *Nat Rev Immunol* 13: 73-74
- Larimore K., McCormick M. W., Robins H. S. and Greenberg P. D. (2012) Shaping of human germline IgH repertoires revealed by deep sequencing. *J Immunol* 189: 3221-3230
- Le S. Q. and Durbin R. (2011) SNP detection and genotyping from low-coverage sequencing data on multiple diploid samples. *Genome Res* 21: 952-960
- Leong J. W., Sullivan R. P. and Fehniger T. A. (2012) Natural Killer Cell Regulation by MicroRNAs in Health and Disease. *J Biomed Biotechnol* 2012: 632329
- Lessard C. J., Li H., Adrianto I., Ice J. A., Rasmussen A., et al. (2013) Variants at multiple loci implicated in both innate and adaptive immune responses are associated with Sjögren's syndrome. *Nat Genet* 45: 1284-1292
- Leung J. M., Davenport M., Wolff M. J., Wiens K. E., Abidi W. M., et al. (2014) IL-22-producing CD4+ cells are depleted in actively inflamed colitis tissue. *Mucosal Immunol* 7: 124-133
- Levin A. M., Iannuzzi M. C., Montgomery C. G., Trudeau S., Datta I., et al. (2013) Association of ANXA11 genetic variation with sarcoidosis in African Americans and European Americans. *Genes Immun* 14: 13-18
- Levy S., Sutton G., Ng P. C., Feuk L., Halpern A. L., et al. (2007) The diploid genome sequence of an individual human. *PLoS Biol* 5: e254
- Lewis B. P., Burge C. B. and Bartel D. P. (2005) Conserved seed pairing, often flanked by adenosines, indicates that candidates of human genes are microRNA targets. *Cell* 120: 15-20
- Li Y., Willer C. J., Ding J., Scheet P. and Abecasis G. R. (2010) MaCH: using sequence and genotype data to estimate haplotypes and unobserved genotypes. *Genet Epidemiol* 34: 816-834
- Li Y., Sidore C., Kang H. M., Boehnke M. and Abecasis G. R. (2011) Low-coverage sequencing: Implications for design of complex trait association studies. *Genome Res* 21: 940-951

- Li Y., Cheng H., Zuo X. B., Sheng Y. J., Zhou F. S., et al. (2013) Association analyses identifying two common susceptibility loci shared by psoriasis and systemic lupus erythematosus in the Chinese Han population. *J Med Genet* 50: 812-818
- Li Z., Chao T. C., Chang K. Y., Lin N., Patil V. S., et al. (2014) The long noncoding RNA THRIL regulates TNFa expression through its interaction with hnRNPL. *Proc Natl Acad Sci U S A* 111: 1002-1007
- Lieberman-Aiden E., van Berkum N. L., Williams L., Imakaev M., Ragozy T., et al. (2009) Comprehensive mapping of long-range interactions reveals folding principles of the human genome. *Science* 326: 289-293
- Lipkin W. I. and Firth C. (2013) Viral surveillance and discovery. *Curr Opin Virol* 3: 199-204
- Litman G. W., Rast J. P. and Fugmann S. D. (2010) The origins of vertebrate adaptive immunity. *Nat Rev Immunol* 10: 543-553
- Liu X., Wang Y., Sun Q., Yan J., Huang J., et al. (2012) Identification of microRNA transcriptome involved in human natural killer cell activation. *Immunol Lett* 143: 208-217
- Loo Y. M. and Gale M., Jr. (2011) Immune signaling by RIG-I-like receptors. *Immunity* 34: 680-692
- Lund F. E. and Randall T. D. (2010) Effector and regulatory B cells: modulators of CD4+ T cell immunity. *Nat Rev Immunol* 10: 236-247
- Lutay N., Ambite I., Gronberg Hernandez J., Rydstrom G., Ragnarsdottir B., et al. (2013) Bacterial control of host gene expression through RNA polymerase II. *J Clin Invest* 123: 2366-2379
- Maecher H. T., Lindstrom T. M., Robinson W. H., Utz P. J., Hale M., et al. (2012) New tools for classification and monitoring of autoimmune diseases. *Nat Rev Rheumatol* 8: 317-328
- Maekawa T., Kufer T. A. and Schulze-Lefert P. (2011) NLR functions in plant and animal immune systems: so far and yet so close. *Nat Immunol* 12: 817-826
- Malinge S., Thiollier C., Chlon T. M., Doré L. C., Diebold L., et al. (2013) Ikaros inhibits megakaryopoiesis through functional interaction with GATA-1 and NOTCH signaling. *Blood* 121: 2440-2451
- Mallal S., Nolan D., Witt C., Masel G., Martin A. M., et al. (2002) Association between presence of HLA-B\*5701, HLA-DR7, and HLA-DQ3 and hypersensitivity to HIV-1 reverse-transcriptase inhibitor abacavir. *Lancet* 359: 727-732
- Marchini J., Howie B., Myers S., McVean G. and Donnelly P. (2007) A new multipoint method for genome-wide association studies by imputation of genotypes. *Nat Genet* 39: 906-913
- Markle J. G., Frank D. N., Mortin-Toth S., Robertson C. E., Feazel L. M., et al. (2013) Sex differences in the gut microbiome drive hormone-dependent regulation of autoimmunity. *Science* 339: 1084-1088
- Marks C. (1983) Immunobiological determinants in organ transplantation. *Ann R Coll Surg Engl* 65: 139-144
- Marks J. D., Hoogenboom H. R., Bonnert T. P., McCafferty J., Griffiths A. D., et al. (1991) Bypassing immunization. Human antibodies from V-gene libraries displayed on phage. *J Mol Biol* 222: 581-597
- Martin J. E., Assassi S., Diaz-Gallo L. M., Broen J. C., Simeon C. P., et al. (2013) A systemic sclerosis and systemic lupus erythematosus pan-meta-GWAS reveals new shared susceptibility loci. *Hum Mol Genet* 22: 4021-4029
- Martinez N. M. and Lynch K. W. (2013) Control of alternative splicing in immune responses: many regulators, many predictions, much still to learn. *Immunol Rev* 253: 216-236
- Maslanik T., Tannura K., Mahaffey L., Loughridge A. B., Bennisson L., et al. (2012) Commensal bacteria and MAMPs are necessary for stress-induced increases in IL-1beta and IL-18 but not IL-6, IL-10 or MCP-1. *PLoS ONE* 7: e50636
- Mavrommatis B., Young G. R. and Kassiotis G. (2013) Counterpoise between the microbiome, host immune activation and pathology. *Curr Opin Immunol* 25: 456-462
- McDonald-McGinn D. M., Fahiminiya S., Revil T., Nowakowska B. A., Suhl J., et al. (2013) Hemizygous mutations in SNAP29 unmask autosomal recessive conditions and contribute to atypical findings in patients with 22q11.2DS. *J Med Genet* 50: 80-90
- McGuckin M. A., Eri R., Simms L. A., Florin T. H. and Radford-Smith G. (2009) Intestinal barrier dysfunction in inflammatory bowel diseases. *Inflamm Bowel Dis* 15: 100-113
- Medvedovic J., Ebert A., Tagoh H., Tamir I. M., Schwickert T. A., et al. (2013) Flexible long-range loops in the VH gene region of the IgH locus facilitate the generation of a diverse antibody repertoire. *Immunity* 39: 229-244
- Meier J., Roberts C., Avent K., Hazlett A., Berrie J., et al. (2013) Fractal organization of the human T cell repertoire in health and after stem cell transplantation. *Biol Blood Marrow Transplant* 19: 366-377
- Meyer E. H., Hsu A. R., Liliental J., Ehrlich A., Florek M., et al. (2013) A distinct evolution of the T-cell repertoire categorizes treatment refractory gastrointestinal acute graft-versus-host disease. *Blood* 121: 4955-4962
- Mikacenic C., Reiner A. P., Holden T. D., Nickerson D. A. and Wurfel M. M. (2013) Variation in the TLR10/TLR1/TLR6 locus is the major genetic determinant of interindividual difference in TLR1/2-mediated responses. *Genes Immun* 14: 52-57
- Miles J. J., Douek D. C. and Price D. A. (2011) Bias in the alpha-beta T-cell repertoire: implications for disease pathogenesis and vaccination. *Immunol Cell Biol* 89: 375-387
- Minton K. (2014) Antiviral immunity: Editing HLA-E expression. *Nat Rev Immunol* 14: 210-211
- Mitelman F., Johansson B. and Mertens F. (2007) The impact of translocations and gene fusions on cancer causation. *Nat Rev Cancer* 7: 233-245
- Montavon T., Soshnikova N., Mascrez B., Joye E., Thevenet L., et al. (2011) A regulatory archipelago controls Hox genes transcription in digits. *Cell* 147: 1132-1145
- Nakano K., Whitaker J. W., Boyle D. L., Wang W. and Firestein G. S. (2013) DNA methylome signature in rheumatoid arthritis. *Ann Rheum Dis* 72: 110-117
- Natoli G. (2010) Maintaining cell identity through global control of genomic organization. *Immunity* 33: 12-24
- Neller M. A., Burrows J. M., Rist M. J., Miles J. J. and Burrows S. R. (2013) High frequency of herpesvirus-specific clonotypes in the human T cell repertoire can remain stable over decades with minimal turnover. *J Virol* 87: 697-700
- Nguyen P., Ma J., Pei D., Obert C., Cheng C., et al. (2011) Identification of errors introduced during high throughput sequencing of the T cell receptor repertoire. *BMC Genomics* 12: 106
- Nishikawa H. and Sakaguchi S. (2014) Regulatory T cells in cancer immunotherapy. *Curr Opin Immunol* 27C: 1-7
- Nolte 't Hoen E. N., Buermans H. P., Waasdorp M., Stoorvogel W., Wauben M. H., et al. (2012) Deep sequencing of RNA from immune cell-derived vesicles uncovers the selective incorporation of small non-coding RNA biotypes with potential regulatory functions. *Nucleic Acids Res* 40: 9272-9285
- Noordermeer D., Branco M. R., Splinter E., Klous P., van Ijcken W., et al. (2008) Transcription and chromatin organization of a housekeeping gene cluster containing an integrated beta-globin locus control region. *PLoS Genet* 4: e1000016
- Northrup D. L. and Zhao K. (2011) Application of ChIP-Seq and related techniques to the study of immune function. *Immunity* 34: 830-842
- Nussenzweig A. and Nussenzweig M. C. (2010) Origin of chromosomal translocations in lymphoid cancer. *Cell* 141: 27-38

- O'Connell R. M., Rao D. S., Chaudhuri A. A. and Baltimore D. (2010) Physiological and pathological roles for microRNAs in the immune system. *Nat Rev Immunol* 10: 111-122
- O'Connor K. S., Parnell G., Patrick E., Ahlenstiel G., Suppiah V., et al. (2014) Hepatic metallothionein expression in chronic hepatitis C virus infection is IFNL3 genotype-dependent. *Genes Immun* 15: 88-94
- Oberg A. L., Kennedy R. B., Li P., Ovsyannikova I. G. and Poland G. A. (2011) Systems biology approaches to new vaccine development. *Curr Opin Immunol* 23: 436-443
- Orrù V., Steri M., Sole G., Sidore C., Viridis F., et al. (2013) Genetic variants regulating immune cell levels in health and disease. *Cell* 155: 242-256
- Oto T., Levvey B. J. and Snell G. I. (2007) Potential refinements of the International Society for Heart and Lung Transplantation primary graft dysfunction grading system. *J Heart Lung Transplant* 26: 431-436
- Ottinger H. D., Ferencik S., Beelen D. W., Lindemann M., Peceny R., et al. (2003) Hematopoietic stem cell transplantation: contrasting the outcome of transplantations from HLA-identical siblings, partially HLA-mismatched related donors, and HLA-matched unrelated donors. *Blood* 102: 1131-1137
- Pagani M., Rossetti G., Panzeri I., de Candia P., Bonnal R. J., et al. (2013) Role of microRNAs and long-non-coding RNAs in CD4(+) T-cell differentiation. *Immunol Rev* 253: 82-96
- Paquette S. G., Banner D., Chi I., Le?n A. J., Xu L., et al. (2014) Pandemic H1N1 influenza A directly induces a robust and acute inflammatory gene signature in primary human bronchial epithelial cells downstream of membrane fusion. *Virology* 448: 91-103
- Peng Z., Cheng Y., Tan B. C., Kang L., Tian Z., et al. (2012) Comprehensive analysis of RNA-Seq data reveals extensive RNA editing in a human transcriptome. *Nat Biotechnol* 30: 253-260
- Petrie H. T. and Zuniga-Pflucker J. C. (2007) Zoned out: functional mapping of stromal signaling microenvironments in the thymus. *Annu Rev Immunol* 25: 649-679
- Pinto S., Michel C., Schmidt-Glenewinkel H., Harder N., Rohr K., et al. (2013) Overlapping gene coexpression patterns in human medullary thymic epithelial cells generate self-antigen diversity. *Proc Natl Acad Sci U S A* 110: E3497-3505
- Pociask D. A., Scheller E. V., Mandalapu S., McHugh K. J., Enelow R. I., et al. (2013) IL-22 Is Essential for Lung Epithelial Repair following Influenza Infection. *Am J Pathol* 182: 1286-1296
- Poland G. A., Ovsyannikova I. G., Jacobson R. M. and Smith D. I. (2007) Heterogeneity in vaccine immune response: the role of immunogenetics and the emerging field of vaccinomics. *Clin Pharmacol Ther* 82: 653-664
- Poland G. A., Ovsyannikova I. G. and Jacobson R. M. (2008) Personalized vaccines: the emerging field of vaccinomics. *Expert Opin Biol Ther* 8: 1659-1667
- Poland G. A., Jacobson R. M. and Ovsyannikova I. G. (2009) Trends affecting the future of vaccine development and delivery: the role of demographics, regulatory science, the anti-vaccine movement, and vaccinomics. *Vaccine* 27: 3240-3244
- Poland G. A., Kennedy R. B. and Ovsyannikova I. G. (2011) Vaccinomics and personalized vaccinology: is science leading us toward a new path of directed vaccine development and discovery? *PLoS Pathog* 7: e1002344
- Poland G. A., Ovsyannikova I. G., Kennedy R. B., Haralambieva I. H. and Jacobson R. M. (2011) Vaccinomics and a new paradigm for the development of preventive vaccines against viral infections. *OMICS* 15: 625-636
- Prokunina-Olsson L., Muchmore B., Tang W., Pfeiffer R. M., Park H., et al. (2013) A variant upstream of IFNL3 (IL28B) creating a new interferon gene IFNL4 is associated with impaired clearance of hepatitis C virus. *Nat Genet* 45: 164-171
- Putintseva E. V., Britanova O. V., Staroverov D. B., Merzlyak E. M., Turchaninova M. A., et al. (2013) Mother and child T cell receptor repertoires: deep profiling study. *Front Immunol* 4: 463
- Quintana-Murci L. and Clark A. G. (2013) Population genetic tools for dissecting innate immunity in humans. *Nat Rev Immunol* 13: 280-293
- Riolobos L., Hirata R. K., Turtle C. J., Wang P. R., Gornalusse G. G., et al. (2013) HLA Engineering of Human Pluripotent Stem Cells. *Mol Ther* 21: 1232-1241
- Robbiani D. F. and Nussenzweig M. C. (2013) Chromosome translocation, B cell lymphoma, and activation-induced cytidine deaminase. *Annu Rev Pathol* 8: 79-103
- Robins H. (2013) Immunosequencing: applications of immune repertoire deep sequencing. *Curr Opin Immunol* 25: 646-652
- Robins H. S., Campregher P. V., Srivastava S. K., Wacher A., Turtle C. J., et al. (2009) Comprehensive assessment of T-cell receptor beta-chain diversity in alphabeta T cells. *Blood* 114: 4099-4107
- Robins H. S., Srivastava S. K., Campregher P. V., Turtle C. J., Andriesen J., et al. (2010) Overlap and effective size of the human CD8+ T cell receptor repertoire. *Sci Transl Med* 2: 47ra64
- Rocha P. P., Micsinai M., Kim J. R., Hewitt S. L., Souza P. P., et al. (2012) Close proximity to Igh is a contributing factor to AID-mediated translocations. *Mol Cell* 47: 873-885
- Rothenberg E. V., Moore J. E. and Yui M. A. (2008) Launching the T-cell-lineage developmental programme. *Nat Rev Immunol* 8: 9-21
- Rothenberg E. V. (2012) Transcriptional drivers of the T-cell lineage program. *Curr Opin Immunol* 24: 132-138
- Roychoudhuri R., Hirahara K., Mousavi K., Clever D., Klebanoff C. A., et al. (2013) BACH2 represses effector programs to stabilize T(reg)-mediated immune homeostasis. *Nature* 498: 506-510
- Sabeti P. C., Reich D. E., Higgins J. M., Levine H. Z., Richter D. J., et al. (2002) Detecting recent positive selection in the human genome from haplotype structure. *Nature* 419: 832-837
- Sajan S. A. and Hawkins R. D. (2012) Methods for identifying higher-order chromatin structure. *Annu Rev Genomics Hum Genet* 13: 59-82
- Sancho D. and Reis e Sousa C. (2012) Signaling by myeloid C-type lectin receptors in immunity and homeostasis. *Annu Rev Immunol* 30: 491-529
- Sandler N. G., Koh C., Roque A., Eccleston J. L., Siegel R. B., et al. (2011) Host response to translocated microbial products predicts outcomes of patients with HBV or HCV infection. *Gastroenterology* 141: 1220-1230, 1230 e1221-1223
- Scheet P. and Stephens M. (2008) Linkage disequilibrium-based quality control for large-scale genetic studies. *PLoS Genet* 4: e1000147
- Schoggins J. W., MacDuff D. A., Imanaka N., Gainey M. D., Shrestha B., et al. (2014) Pan-viral specificity of IFN-induced genes reveals new roles for cGAS in innate immunity. *Nature* 505: 691-695
- Schroder K. and Tschopp J. (2010) The inflammasomes. *Cell* 140: 821-832
- Schwab C., Berry D., Rauch I., Rennisch I., Ramesmayer J., et al. (2014) Longitudinal study of murine microbiota activity and interactions with the host during acute inflammation and recovery. *ISME J*
- Seed K. D., Lazinski D. W., Calderwood S. B. and Camilli A. (2013) A bacteriophage encodes its own CRISPR/Cas adaptive response to evade host innate immunity. *Nature* 494: 489-491
- Shalek A. K., Satija R., Adiconis X., Gertner R. S., Gaublomme J. T., et al. (2013) Single-cell transcriptomics reveals bimodality in expression and splicing in immune cells. *Nature* 498: 236-240
- Shay T. and Kang J. (2013) Immunological Genome Project and systems immunology. *Trends Immunol* 34: 602-609

- Sherwood A. M., Emerson R. O., Scherer D., Habermann N., Buck K., et al. (2013) Tumor-infiltrating lymphocytes in colorectal tumors display a diversity of T cell receptor sequences that differ from the T cells in adjacent mucosal tissue. *Cancer Immunol Immunother* 62: 1453-1461
- Snyder T. M., Khush K. K., Valantine H. A. and Quake S. R. (2011) Universal noninvasive detection of solid organ transplant rejection. *Proceedings of the National Academy of Sciences of the United States of America* 108: 6229-6234
- Starzl R., Brandacher G., Lee W. P., Carbonell J., Zhang W., et al. (2013) Review of the early diagnoses and assessment of rejection in vascularized composite allotransplantation. *Clin Dev Immunol* 2013: 402980
- Steinert E. M., Schwartz R. H. and Singh N. J. (2012) At low precursor frequencies, the T-cell response to chronic self-antigen results in anergy without deletion. *Eur J Immunol* 42: 2875-2880
- Su X., Qian C., Zhang Q., Hou J., Gu Y., et al. (2013) miRNomes of haematopoietic stem cells and dendritic cells identify miR-30b as a regulator of Notch1. *Nat Commun* 4: 2903
- Sui J., Hwang W. C., Perez S., Wei G., Aird D., et al. (2009) Structural and functional bases for broad-spectrum neutralization of avian and human influenza A viruses. *Nat Struct Mol Biol* 16: 265-273
- Swaminathan S., Hu X., Zheng X., Kriga Y., Shetty J., et al. (2013) Interleukin-27 treated human macrophages induce the expression of novel microRNAs which may mediate anti-viral properties. *Biochem Biophys Res Commun* 434: 228-234
- Takeuchi O. and Akira S. (2010) Pattern recognition receptors and inflammation. *Cell* 140: 805-820
- Tamouza R., El Kassar N., Schaeffer V., Carbonnelle E., Tatarski Z., et al. (2000) A novel HLA-B\*39 allele (HLA-B\*3916) due to a rare mutation causing cryptic splice site activation. *Hum Immunol* 61: 467-473
- Tewhey R., Bansal V., Torkamani A., Topol E. J. and Schork N. J. (2011) The importance of phase information for human genomics. *Nat Rev Genet* 12: 215-223
- Thomas R., Apps R., Qi Y., Gao X., Male V., et al. (2009) HLA-C cell surface expression and control of HIV/AIDS correlate with a variant upstream of HLA-C. *Nat Genet* 41: 1290-1294
- Tishkoff S. A., Dietzsch E., Speed W., Pakstis A. J., Kidd J. R., et al. (1996) Global patterns of linkage disequilibrium at the CD4 locus and modern human origins. *Science* 271: 1380-1387
- Tokoyoda K., Hauser A. E., Nakayama T. and Radbruch A. (2010) Organization of immunological memory by bone marrow stroma. *Nat Rev Immunol* 10: 193-200
- Trautmann L., Rimbart M., Echasserieau K., Saulquin X., Neveu B., et al. (2005) Selection of T cell clones expressing high-affinity public TCRs within Human cytomegalovirus-specific CD8 T cell responses. *J Immunol* 175: 6123-6132
- Tsai A. G., Lu H., Raghavan S. C., Muschen M., Hsieh C. L., et al. (2008) Human chromosomal translocations at CpG sites and a theoretical basis for their lineage and stage specificity. *Cell* 135: 1130-1142
- Tsoi L. C., Spain S. L., Knight J., Ellinghaus E., Stuart P. E., et al. (2012) Identification of 15 new psoriasis susceptibility loci highlights the role of innate immunity. *Nat Genet* 44: 1341-1348
- Turchaninova M. A., Britanova O. V., Bolotin D. A., Shugay M., Putintseva E. V., et al. (2013) Pairing of T-cell receptor chains via emulsion PCR. *Eur J Immunol* 43: 2507-2515
- Vahedi G., Takahashi H., Nakayama S., Sun H. W., Sartorelli V., et al. (2012) STATs Shape the Active Enhancer Landscape of T Cell Populations. *Cell* 151: 981-993
- Venturi V., Price D. A., Douek D. C. and Davenport M. P. (2008) The molecular basis for public T-cell responses? *Nat Rev Immunol* 8: 231-238
- Viatte S., Plant D. and Raychaudhuri S. (2013) Genetics and epigenetics of rheumatoid arthritis. *Nat Rev Rheumatol* 9: 141-153
- Vigano M. A., Ivanek R., Balwierz P., Berninger P., van Nimwegen E., et al. (2013) An epigenetic profile of early T-cell development from multipotent progenitors to committed T-cell descendants. *Eur J Immunol*
- Vladimer G. I., Marty-Roix R., Ghosh S., Weng D. and Lien E. (2013) Inflammasomes and host defenses against bacterial infections. *Curr Opin Microbiol* 16: 23-31
- Walsh D., McCarthy J., O'Driscoll C. and Melgar S. (2013) Pattern recognition receptors--molecular orchestrators of inflammation in inflammatory bowel disease. *Cytokine Growth Factor Rev* 24: 91-104
- Wang C., Sanders C. M., Yang Q., Schroeder H. W., Jr., Wang E., et al. (2010) High throughput sequencing reveals a complex pattern of dynamic interrelationships among human T cell subsets. *Proc Natl Acad Sci U S A* 107: 1518-1523
- Wang C., Krishnakumar S., Wilhelmy J., Babrzadeh F., Stepanyan L., et al. (2012) High-throughput, high-fidelity HLA genotyping with deep sequencing. *Proc Natl Acad Sci U S A* 109: 8676-8681
- Wang C., Ahlford A., Laxman N., Nordmark G., Eloranta M. L., et al. (2013) Contribution of IKBKE and IFIH1 gene variants to SLE susceptibility. *Genes Immun* 14: 217-222
- Wang P., Gu Y., Zhang Q., Han Y., Hou J., et al. (2012) Identification of resting and type I IFN-activated human NK cell miRNomes reveals microRNA-378 and microRNA-30e as negative regulators of NK cell cytotoxicity. *J Immunol* 189: 211-221
- Wang X., Lin Z., Gao L., Wang A., Wan Z., et al. (2013) Exome sequencing reveals a signal transducer and activator of transcription 1 (STAT1) mutation in a child with recalcitrant cutaneous fusariosis. *J Allergy Clin Immunol* 131: 1242-1243
- Wang X., Wang H. K., Li Y., Hafner M., Banerjee N. S., et al. (2014) microRNAs are biomarkers of oncogenic human papillomavirus infections. *Proc Natl Acad Sci U S A* 111: 4262-4267
- Warren E. H., Matsen F. A. t. and Chou J. (2013) High-throughput sequencing of B- and T-lymphocyte antigen receptors in hematology. *Blood* 122: 19-22
- Warren R. L., Freeman J. D., Zeng T., Choe G., Munro S., et al. (2011) Exhaustive T-cell repertoire sequencing of human peripheral blood samples reveals signatures of antigen selection and a directly measured repertoire size of at least 1 million clonotypes. *Genome Res* 21: 790-797
- Weaver W. M., Tseng P., Kunze A., Masaeli M., Chung A. J., et al. (2014) Advances in high-throughput single-cell microtechnologies. *Curr Opin Biotechnol* 25C: 114-123
- Wei G., Wei L., Zhu J., Zang C., Hu-Li J., et al. (2009) Global mapping of H3K4me3 and H3K27me3 reveals specificity and plasticity in lineage fate determination of differentiating CD4+ T cells. *Immunity* 30: 155-167
- Wei Y., Chen R., Dimicoli S., Bueso-Ramos C., Neuberger D., et al. (2013) Global H3K4me3 genome mapping reveals alterations of innate immunity signaling and overexpression of JMJD3 in human myelodysplastic syndrome CD34+ cells. *Leukemia* 27: 2177-2186
- Wherry E. J. (2011) T cell exhaustion. *Nat Immunol* 12: 492-499
- Whisnant A. W., Bogerd H. P., Flores O., Ho P., Powers J. G., et al. (2013) In-depth analysis of the interaction of HIV-1 with cellular microRNA biogenesis and effector mechanisms. *MBio* 4: e000193
- Wie S. H., Du P., Luong T. Q., Rought S. E., Beliakova-Bethell N., et al. (2013) HIV Downregulates Interferon-Stimulated Genes in Primary Macrophages. *J Interferon Cytokine Res* 33: 90-95
- Wilson P. C. and Andrews S. F. (2012) Tools to therapeutically harness the human antibody response. *Nat Rev Immunol* 12: 709-719

- Winters G. L. and McManus B. M. (1996) Consistencies and controversies in the application of the International Society for Heart and Lung Transplantation working formulation for heart transplant biopsy specimens. Rapamycin Cardiac Rejection Treatment Trial Pathologists. *J Heart Lung Transplant* 15: 728-735
- Woodsworth D. J., Castellarin M. and Holt R. A. (2013) Sequence analysis of T-cell repertoires in health and disease. *Genome Med* 5: 98
- Wu D., Sherwood A., Fromm J. R., Winter S. S., Dunsmore K. P., et al. (2012) High-throughput sequencing detects minimal residual disease in acute T lymphoblastic leukemia. *Sci Transl Med* 4: 134ra163
- Wu X., Zhou T., Zhu J., Zhang B., Georgiev I., et al. (2011) Focused evolution of HIV-1 neutralizing antibodies revealed by structures and deep sequencing. *Science* 333: 1593-1602
- Yamane A., Robbiani D. F., Resch W., Bothmer A., Nakahashi H., et al. (2013) RPA Accumulation during Class Switch Recombination Represents 5'-3' DNA-End Resection during the S-G2/M Phase of the Cell Cycle. *Cell Rep* 3: 138-147
- Yang G., Yang L., Zhao Z., Wang J. and Zhang X. (2012) Signature miRNAs Involved in the Innate Immunity of Invertebrates. *PLoS ONE* 7: e39015
- Yu Z., Garner C., Zlogas A., Anton-Culver H. and Schaid D. J. (2009) Genotype determination for polymorphisms in linkage disequilibrium. *BMC Bioinformatics* 10: 63
- Yurkovetskiy L., Burrows M., Khan A. A., Graham L., Volchkov P., et al. (2013) Gender bias in autoimmunity is influenced by microbiota. *Immunity* 39: 400-412
- Zhang J. A., Mortazavi A., Williams B. A., Wold B. J. and Rothenberg E. V. (2012) Dynamic Transformations of Genome-wide Epigenetic Marking and Transcriptional Control Establish T Cell Identity. *Cell* 149: 467-482
- Zhang Y., Gostissa M., Hildebrand D. G., Becker M. S., Boboila C., et al. (2010) The role of mechanistic factors in promoting chromosomal translocations found in lymphoid and other cancers. *Adv Immunol* 106: 93-133
- Zhao Z., Tavoosidana G., Sjolinder M., Gondor A., Mariano P., et al. (2006) Circular chromosome conformation capture (4C) uncovers extensive networks of epigenetically regulated intra- and interchromosomal interactions. *Nat Genet* 38: 1341-1347
- Zheng X., Shen J., Cox C., Wakefield J. C., Ehm M. G., et al. (2014) HIBAG-HLA genotype imputation with attribute bagging. *Pharmacogenomics J* 14: 192-200
- Zheng Y., Zha Y., Spaapen R. M., Mathew R., Barr K., et al. (2013) Egr2-dependent gene expression profiling and ChIP-Seq reveal novel biologic targets in T cell anergy. *Mol Immunol* 55: 283-291
- Zhou L., Chong M. M. and Littman D. R. (2009) Plasticity of CD4+ T cell lineage differentiation. *Immunity* 30: 646-655



FOR RESEARCH USE ONLY

© 2014 Illumina, Inc. All rights reserved. Illumina, IlluminaDx, BaseSpace, BeadArray, BeadXpress, cBot, CSPRO, DASL, DesignStudio, Eco, GAllx, Genetic Energy, Genome Analyzer, GenomeStudio, GoldenGate, HiScan, HiSeq, Infinium, iSelect, MiSeq, Nextera, NuPCR, SeqMonitor, Solexa, TruSeq, TruSight, VeraCode, the pumpkin orange color, and the Genetic Energy streaming bases design are trademarks or registered trademarks of Illumina, Inc. All other brands and names contained herein are the property of their respective owners.  
Pub. No. 470-2014-003 Current as of 24 July 2014

